

PICARRO

G2207-*i* Analyzer for O₂ and $\delta^{18}\text{O}$ in Air User Manual



Picarro Inc.

3105 Patrick Henry Drive
Santa Clara, CA 95054, USA

Phone: +1 408 962 3944

www.picarro.com

Document Number 40-0034 Rev. B

Picarro Analyzer User Manual

Thank you for purchasing a Picarro product. Your Picarro system is a quality product that has been designed and manufactured to provide reliable performance.

This User Manual (UM) is an important part of your purchase as it will help familiarize you with the system and explain the numerous features that have been designed into it. Please read this manual thoroughly before using your Picarro system.

Please contact Picarro or your authorized Picarro distributor should you have questions regarding specific applications or if you require additional information.

NOTICES

DISCLAIMER AND RESERVATION OF RIGHTS

Picarro has prepared this manual solely for the information and use by its customers as a guide for the selection, installation, operation, and maintenance of the products described.

EXCEPT AS PROVIDED IN THE TERMS AND CONDITIONS OF SALE FOR PICARRO PRODUCTS, PICARRO ASSUMES NO LIABILITY WHATSOEVER, AND PICARRO DISCLAIMS ANY EXPRESS OR IMPLIED WARRANTY, RELATING TO SALE OR USE OF PICARRO PRODUCTS, INCLUDING LIABILITY OR WARRANTIES RELATING TO MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.

Picarro reserves the right to change or update the contents of this manual and the specifications of its products at any time, without notice. Picarro has endeavored to include information that is current and accurate as of the date of the publication or revision of this document, but Picarro does not guarantee that this document is error free or that it is accurate with regard to any particular specification.

Picarro expressly reserves all intellectual property rights, including all intellectual property rights relating to any product described in this manual. This document does not grant any license, express or implied, by estoppel or otherwise, to any intellectual property rights of Picarro or any third party.

PATENTS

The products described in this manual are subject to Picarro patents and patents pending. Information about Picarro patents applicable to these products is available at www.picarro.com/company/patents

TRADEMARKS

Picarro and the Picarro logo are trademarks of Picarro, Inc.

Windows® is a trademark of Microsoft Corporation

Swagelok® is a trademark of Swagelok Company

Copyright © 2020 Picarro, Inc. All rights reserved

Contact Information

General Technical Support:

Email: support@picarro.com

Phone: +1 408 962 3991

European Technical Center:

Email: support@picarro.com

Phone: +31 85 888 1650

Customer Service:

Email: orders@picarro.com

Phone: +1 408 962 3992

© 2020 Picarro, Inc. All rights reserved.
3105 Patrick Henry Drive
Santa Clara, California, 95054 USA
Phone 408.962.3900 • Fax 408.962.3200
www.picarro.com

Table of Contents

1. INTRODUCTION TO TECHNOLOGY	7
1.1 Cavity Ring-Down Spectroscopy (CRDS).....	7
1.2 Relating ring-down time to absorption intensity	8
1.3 Converting absorption intensity to concentration.....	9
1.4 Spectral precision and high sensitivity measurements	10
2. CONVENTIONS	12
3. ACRONYMS.....	13
4. SAFETY	14
4.1 General Safety	14
4.2 Laser Safety	15
5. UNPACKING THE ANALYZER	17
5.1 Inspect the Shipping Boxes	17
5.2 Unpack the Shipping Boxes	17
6. ANALYZER OVERVIEW	20
6.1 Intended Use	20
6.2 External Vacuum Pump.....	20
6.3 Analyzer Specifications	21
7. INSTALLATION	22
7.2 Analyzer and Vacuum Setup	24
7.3 Connecting to the Analyzer Inlet.....	25
7.4 Setting Up a Monitor, Keyboard, and Mouse.....	26
8. ANALYZER BASIC OPERATION.....	28
8.1 Startup Procedure	28
8.2 Desktop Icons and Folders.....	29
8.3 Shutdown Procedure	31
8.4 In case of an Electrical Power Outage	32

9. LIST OF GUI FUNCTIONS.....	34
9.1 GUI Overview	34
9.2 Settings, Tools, and Help Menus.....	34
9.3 Alarm Panel.....	35
9.4 Instrument status.....	36
9.5 Data Log Filename and Path.....	36
9.6 Data Window	37
9.7 Status Log Window	40
10. FILE MANAGEMENT.....	42
10.1 Overview	42
10.2 File Archival and Automatic Deletion of Old Files.....	44
10.3 Set Up Tool	46
11. CALIBRATION.....	53
11.1 Slope and Offset.....	53
11.2 Calibration Methodology.....	53
11.3 Calibration Software Utility tools.....	56
12. TROUBLESHOOTING	65
12.1 Power LED on Analyzer Does Not Illuminate	65
12.2 User Interface Program Does Not Start.....	65
12.3 Sample Pressure Cannot Be Controlled to Appropriate Value for Concentration Measurements	66
12.4 User Interface Program “Freezes”/Does Not Update Graphs as Data are Collected.	66
13. SERVICE AND MAINTENANCE.....	67
13.1 Particulate Filter Replacement	67
13.2 Cleaning	71
14. TRANSPORTATION AND STORAGE.....	72
15. LIMITED WARRANTY	73
APPENDIX A – REMOTE DATA ACCESS	75
A.1 Picarro Serial Communication.....	75

A.2 Remote Data Access	75
APPENDIX B – THE DATA FILE VIEWER	78
B.1 Introduction	78
B.2 The File Menu	78
B.3 New – Time Series Plot	83
B.4 Time Series Viewer Menus	83
APPENDIX C – EXTERNAL VALVE SEQUENCER.....	92
C.1 Introduction	92
C.2 Default Configuration	92
C.3 Setting Up Solenoid Valves	92
C.4 Setting Up Rotary Selector Valve	93
C.5 Valve Sequencer Software	94
C.6 Configuring a Valve Sequence	95
C.7 Loading and Running a Saved Sequence	95

1. INTRODUCTION TO TECHNOLOGY

Picarro analyzers use time-based, optical absorption spectroscopy of the target gases to determine concentration in a sample. They are based on wavelength-scanned cavity ring-down spectroscopy (WS-CRDS), a technology in which light travels many times through the sample, creating a very long effective path length for the light to interact with the target gas, thus enabling excellent detection sensitivity in a compact and rugged instrument.

The Picarro analyzer is comprised of two modules:

- The **Analyzer** contains the spectrometer, sample chamber, and a computer with a hard drive to store and analyze data. The single analyzer module controls the operation of the system and converts spectroscopic measurements into gas concentration data.
- The **External Vacuum Pump** draws the sample gas through the instrument.

1.1 Cavity Ring-Down Spectroscopy (CRDS)

Nearly every small gas-phase molecule (e.g., CO₂, H₂O, H₂S, NH₃) and isotopologue (e.g., H₂¹⁸O, ¹³CO₂, ¹⁵N¹⁴N¹⁶O) uniquely absorb specific wavelengths of near-infrared light. The strength of the light absorption is related to the concentration of a molecule in a sample and the distance that light travels through the sample, called the path length.


Conventional infrared spectrometers are typically only sensitive enough to detect trace gases at levels in the part-per-million. Cavity Ring-Down Spectroscopy (CRDS), on the other hand, is one thousand to one million more times sensitive.

The increased sensitivity of CRDS is due to the design of the sample cavity and the time-based measurement. In the cavity, a series of mirrors reflects the infrared light through the sample, increasing the path length. For a Picarro cavity of only 25 cm in length, the effective path length of the cavity can be over 20 kilometers.

In Picarro analyzers, light from a single-frequency laser enters a cavity where three mirrors reflect the laser light as seen in Figure 1. The light enters through the mirror closest to the laser, bounces off the angled mirror in the lower right corner of the cavity, travels to the hemispherical mirror at the top of the cavity, bounces toward the mirror in the lower left corner of the cavity, and then returns to the first mirror. This motion becomes a

continuous traveling light wave, which is represented by the dark orange path in Figure 1.

Figure 1
Schematic of
the Picarro
CRDS analyzer
cavity

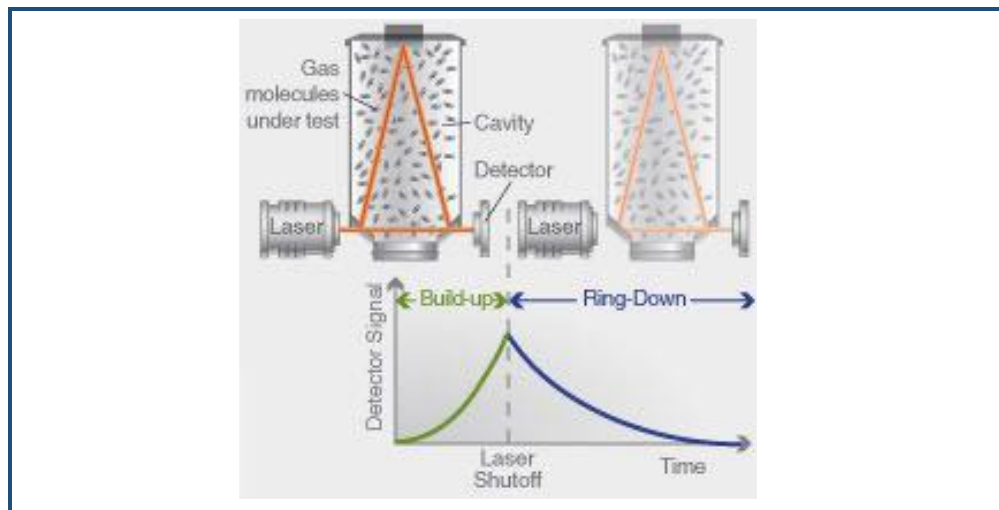

When the laser is on, the cavity quickly fills with laser light. A small amount of the laser light is transmitted through the mirror closest to the photodetector, which turns the incident light into a signal that is directly proportional to the light intensity in the cavity.

When the photodetector signal reaches a threshold level (in a few tens of microseconds), the laser is turned off. The light contained within the cavity continues to bounce between the mirrors (about 40,000 times). Since the mirrors have slightly less than 100% reflectivity (99.999%), the light inside the cavity steadily leaks out of the cavity. The intensity of the light reaching the detector decreases, falling exponentially until it reaches zero. This decay, or "ring-down," is measured in real time by the photodetector.

1.2 Relating ring-down time to absorption intensity


The time it takes to ring-down is inversely related to the total optical loss in the cavity, including the strength of molecular absorption at a given wavelength of light. For an empty cavity, the time it takes for the intensity to decrease by a given percent is determined solely by the reflectivity of the mirrors. A cavity containing gas that absorbs light will have a shorter ring-down time than an empty cavity. As the light circulates in a cavity with a gas sample, the molecular absorption by the gas results in a decrease of the light intensity.

Determining absorption intensity at a specific wavelength requires comparing the ring-down time of an empty cavity to the ring-down time of a cavity that contains gas. A cavity can be empty if it contains no gas; it will

also appear empty if the molecules of the sample inside the cavity do not interact with the specific wavelength of light.

Picarro instruments gather measurements from an “empty” cavity by switching the light to wavelengths that are not absorbed by the target molecules. The analyzer subsequently measures ring-down times at wavelengths that are absorbed by the target gas. The analyzer automatically and continuously compares these two types of ring-down times, and the software uses those comparisons to calculate absorption intensities.

Figure 1:
Light intensity
as a function
of time in a
CRDS system

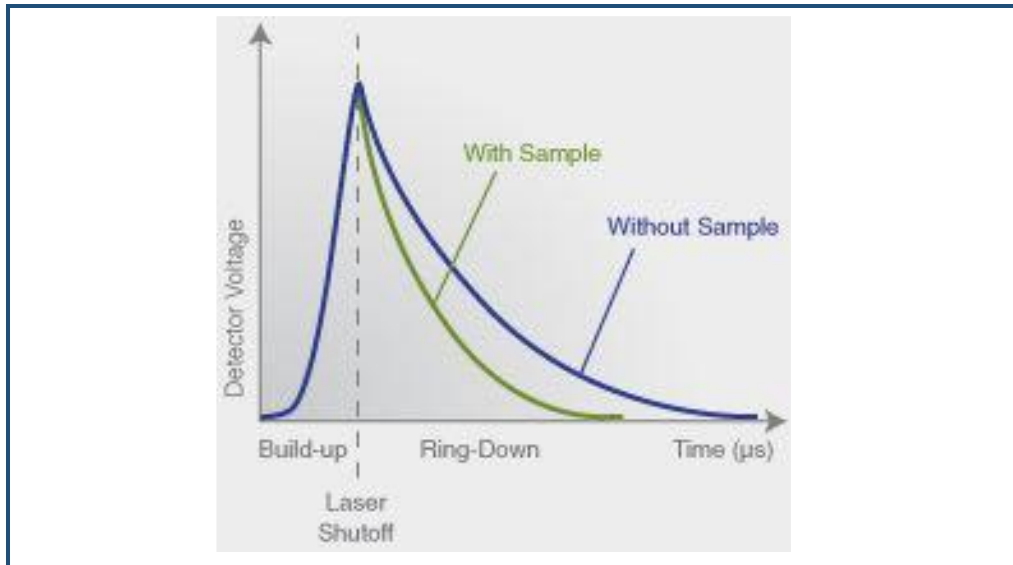


1.3 Converting absorption intensity to concentration

Plotting the absorbance at each measured wavelength generates an optical spectrum. This spectrum contains absorbance peaks that are unique to each molecule in the sample. The height of a particular absorption peak is proportional to the concentration of a molecule that generated the signal.

The height of the peak is calculated by subtracting the maximal absorbance from the baseline absorbance. Figure 3 shows a plot of ideal optical spectra with a clean, uniform baseline on either side of the absorption peak.

Figure 2:
Absorption
Spectral Curve

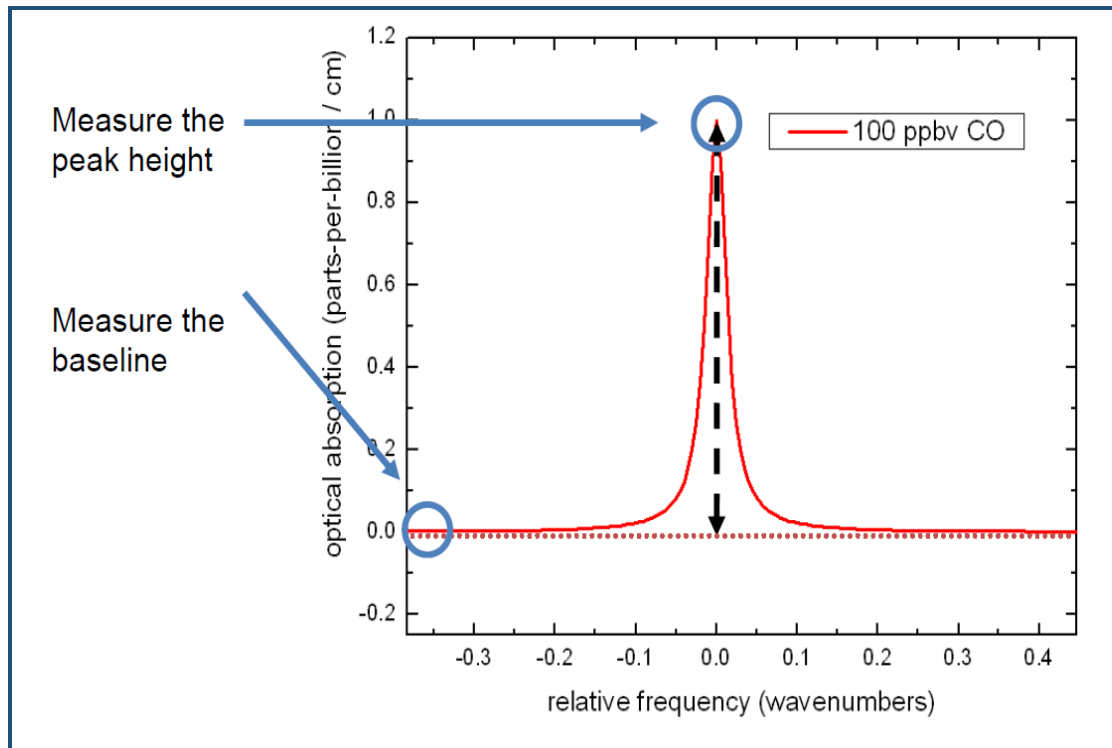

However, optical spectra often contain several absorption lines, nested closely together. A particular absorption peak may be visible between lines, but the absorption may not return to the baseline before it rises in response to another molecule.

Picarro analyzers calculate the baseline underneath a poorly resolved peak by modeling the absorption peaks from other surrounding molecules and subtracting contributions from neighboring peaks to the absorption intensity.

1.4 Spectral precision and high sensitivity measurements

Picarro analyzers contain two features that provide high spectral precision:

- A proprietary **wavelength monitor (WLM)** that measures the absolute laser wavelength to a precision that is a few orders of magnitude narrower than the spectral linewidth: Picarro's patented WLM measures absolute laser wavelength to a precision more than 1,000 times narrower than the observed Doppler-broadened linewidth for small gas-phase molecules. The instruments lock the laser to the WLM, and then the monitor tunes to wavelengths known to be maximally and minimally absorbed by the target molecule. The result is closely clustered absorption intensities, measured at wavelengths just before peak

absorption, at peak absorption, and just after peak absorption, as the absorbance returns to the baseline.







- Precise **temperature and pressure control** in the sample cavity: Accurate absorption measurements at precisely known wavelengths account for little unless the temperature and pressure of the CRDS measurement cavity are known. The observed line intensity and shape depend on the temperature and pressure inside the sample cavity. Small temperature and pressure instabilities can result in large concentration errors due to fluctuating peak heights and baselines. To completely minimize instrument measurement drift, temperature and pressure must be actively stabilized to constant values.

For precise temperature control, the sample cavity is surrounded by layers of thermally insulating material to provide a high degree of passive thermal stability. The cavity is further actively stabilized by means of a solid-state heating system locked to the output of a thermal sensor. This enables the temperature of the cavity to be within 20 mK of the set temperature.

For precise pressure control, the cavity pressure is monitored using a high-linearity pressure transducer. The system computer uses this pressure data in a feedback loop to control proportional valves that adjust the inlet and outlet gas flow of the cavity.

2. CONVENTIONS

The following icons are used throughout this manual to emphasize important information in the text. These icons indicate dangers to either the operator or to the analyzer, and other important information.

	<p>Consult the user's manual for important information (When you see this symbol placed at hazard points on equipment, consult the user manual).</p>
	<p>NOTE is important information that you should be aware of before proceeding.</p>
 <p>WARNING</p>	<p>LASER WARNING alerts you of a laser danger.</p>
 <p>CAUTION</p>	<p>CAUTION alerts you of a potential danger to equipment or to the user.</p>
 <p>WARNING</p>	<p>WARNING indicates an imminent danger to the user.</p>
 <p>REMINDER</p>	<p>REMINDER is a helpful hint to procedures listed in the text.</p>

3. ACRONYMS

This manual includes various acronyms. For definitions, see below:

Acronym	Definition
CRDS	Cavity Ring-Down Spectroscopy
GUI	Graphical User Interface
cm	centimeters
mm	millimeters
" (as in 1/4")	Inches
O ₂	Oxygen
NH ₃	Ammonia
H ₂ O	Water
HB	Hotbox
WB	Warm box
ppm	Parts Per Million
ppb	Parts Per Billion
‰	per mil
°C	degrees Celsius

4. SAFETY

4.1 General Safety

CDRH Certification

This Picarro Analyzer complies with 21 CFR Chapter 1, sub-chapter J, and is classified as a Class 1 laser system when all panels and covers are on.

CE Certification

This Picarro Analyzer complies with the European standards and the instrument is affixed with a CE label. This CE label is located on the rear of the instrument.



WARNING

Using this analyzer in a manner not specified by Picarro may result in damage to the analyzer and render it unsafe to operate.



WARNING

This analyzer is for indoor use only and has an ingress protection rating of IPx-0. Analyzer is NOT protected against exposure to water including dripping, spraying, splashing or immersion.



WARNING

Do not operate in an explosive atmosphere! Do not operate in the presence of flammable gases or fumes.



CAUTION

The analyzer contains no user serviceable components except the particulate filter. Do not attempt repairs; instead, report all problems to Picarro Customer Service or your local distributor. Please contact Picarro if you have any questions regarding the safe operation of this equipment.



CAUTION

The inlet gas connector on the back panel of the Analyzer, and its immediate vicinity, runs hot during operation of the analyzer. Take care when connecting gas lines or working at the rear of the instrument to wear protective gloves or avoid contact with these surfaces.



WARNING

This analyzer weighs 46 lbs. (20.9 kg). Use the technique described below when lifting the analyzer.

- a. Before lifting, inspect the unit for slippery substances or sharp edges.
 - b. Lift with two people, one on each side of the analyzer.
 - c. Crouch down and stay close to the unit. Always keep your back as straight as possible.
 - d. Position your feet for sturdy balance. Lift with your legs, not your back.
 - e. Do not twist the back while carrying the unit. Rotate direction with hip joints.
 - f. Lower the unit by bending at the knees.
-

4.2 Laser Safety



WARNING

This equipment is classified as a Class 1 laser product with an embedded 3B laser in accordance with EN 60825-1:2014. Do not to open the enclosure where this label is placed; there are no user serviceable parts inside.

The following Laser Safety Label is affixed to the outer cover of the analyzer.



WARNING

The laser is a Class3B when exposed. Only operate or service this device in accordance with the instructions in this guide, and only open the device in an approved laser safe service area using appropriate laser-safety glasses.

The following **Laser Safety Label** is affixed to the inside of the analyzer:

Figure 4:
Laser Safety
Label



WARNING

Use of controls or adjustments or performance of procedures other than those specified herein may result in hazardous radiation exposure.

5. UNPACKING THE ANALYZER

5.1 Inspect the Shipping Boxes

Picarro products are inspected and tested before leaving the factory. Their packing containers have been designed to keep the equipment safe from damage during transit.

Inspect the condition of the boxes upon arrival. The larger box includes the analyzer and most of the accessories. Even if the outer box shows damage, the inner box holding the analyzer is cushioned enough that it will protect the instrument under most circumstances.

If the equipment does appear to be damaged, photograph the damages and contact Picarro (email pictures if possible) as soon as possible.

5.2 Unpack the Shipping Boxes


This section describes the contents of the shipping boxes:

- Inspect each item to ensure it is not damaged.
- If items are missing, contact Picarro.
- Keep the shipping packages to reuse when transporting the analyzer.
- Contact Picarro for options on transporting systems to remote labs.



This analyzer weighs 46 lbs. (20.9 kg). Use the technique outlined in the General Safety section on page 15 when lifting or moving the analyzer.

Figure 3:
Box Contents

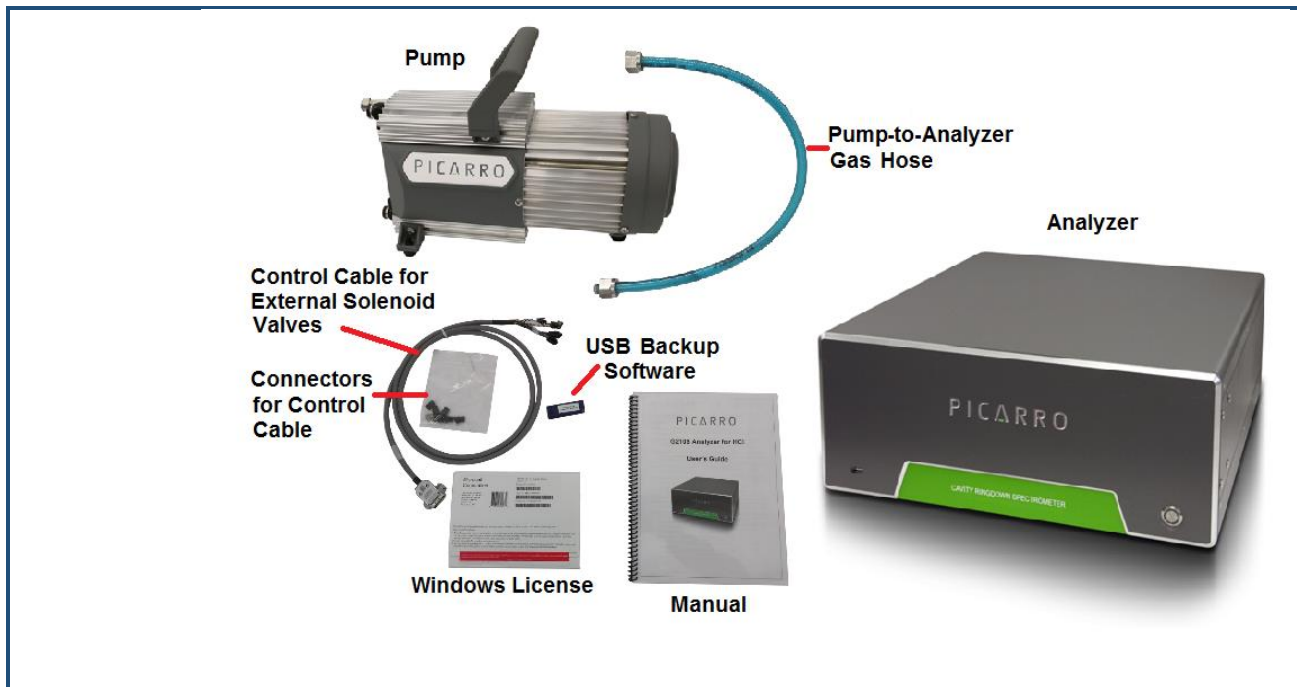

Box One: Analyzer and Accessories

Item (qty)	Description
Analyzer (1)	Includes all the data acquisition, control, and communications hardware and firmware to perform all gas handling, spectral collection and analysis.
A/C Power Cables (1)	A power cable with connectors appropriate to your country is provided. The analyzer automatically adjusts to local voltage.
Keyboard (1)	USB keyboard
Mouse (1)	USB mouse
Control Cable (1)	For External Solenoid Valves
Nut (1) and Ferrules (2)	For connecting input line to analyzer INPUT
Document Packet (1)	Includes this manual, certificate of compliance, and Windows License.
USB Flash Drive	Contains backup software.

Box Two: Vacuum Pump and Accessories

Item (qty)	Description
Pump (1)	Provides vacuum required for sample gas sequencing into and out of the analyzer.
A/C Power Cable (1)	A power cable with connectors appropriate to your country is provided. The pump automatically adjusts to local voltage.
Vacuum Hose (1)	Hose to connect the pump to the analyzer.
Pump Manual (1)	Detailed instructions for pump.

6. ANALYZER OVERVIEW


6.1 Intended Use

The G2207-*i* analyzer measures concentrations of O₂ and δ¹⁸O in air using Picarro's patented Cavity Ring-Down Spectroscopy (CRDS). The analyzer can be deployed for air monitoring applications in a lab or in the field, allowing in-situ analysis of trace and ambient amounts of O₂ and δ¹⁸O.

6.2 External Vacuum Pump

The external vacuum pump is used to maintain cavity pressure inside the analyzer. The pump should be connected and running whenever the analyzer is in use.

Figure 4:
External
vacuum pump

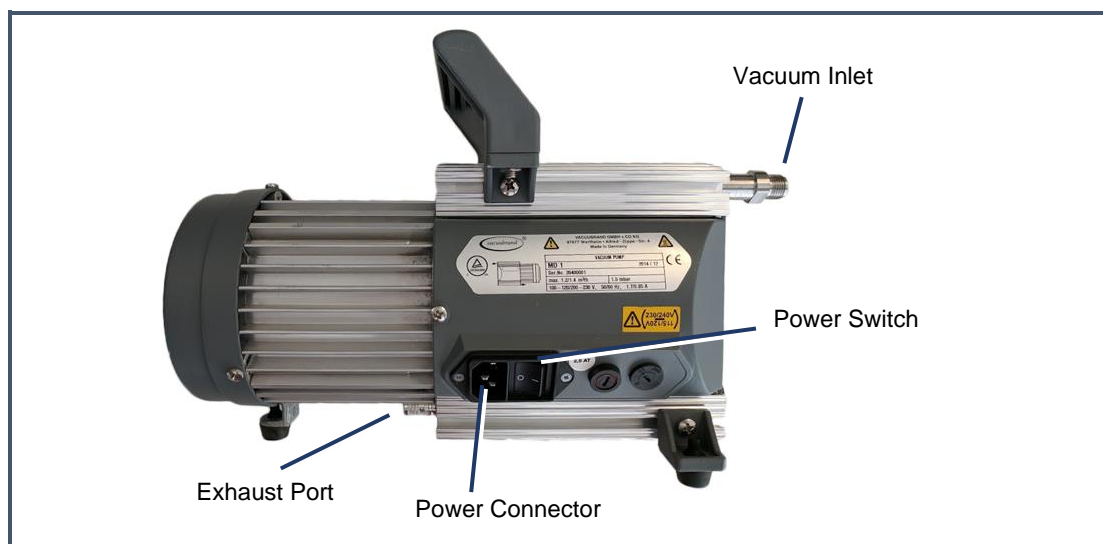

6.3 Analyzer Specifications

Weight (Total)	60.4 lbs. (27.4 kg), Including external pump
Analyzer	46 lbs. (20.9 kg)
Pump	14.4 lbs. (6.5 kg)
Analyzer Dimensions	Length: 17.55" (43.2 cm) Width: 17" (44.6 cm) Height: 7" (17.8 cm)
Temperature Range	Storage: -10 °C to 50 °C; Operation: 10 °C to 35 °C
Sample Flowrate	80 to 110 sccm at 760 torr (101 kPa)
Ambient Humidity Range	< 85% R.H. non-condensing
Maximum Altitude	10,000 ft. (operation)
Clearance	Front: 6" (15.3 cm); Rear: 6" (15.3 cm)
Power Requirements	100 to 240 VAC; 47 to 63 Hz (auto-sensing)
Startup Power	<375 W at start-up (Analyzer and Pump)
Steady-state Power	120 W (Analyzer), 150 W (Pump) Steady-state operation
Mains Supply Voltage Fluctuation	±10% of the nominal voltage
Minimum Rated Circuit Amperage	10 A @115 VAC 5 A @230 VAC
Liquid Ingress Protection	None

7. INSTALLATION

This section describes the setup and installation of the Picarro Analyzer. Please read and understand this section thoroughly before proceeding with the installation.



WARNING

Using this analyzer in a manner not specified by Picarro may result in damage to the analyzer and render it unsafe to operate.



WARNING

Do not attach electrical power to or start the analyzer until *after* attaching and turning on the External Vacuum Pump. Do not disconnect the vacuum line while the analyzer is running. Failure to do so could result in damage to the optics.



WARNING

Picarro sells USB enabled devices, such as GPS, which is approved for use. Please do not connect USB hubs or unapproved USB devices, other than flash drives to the computer because they can interfere with the operation of the analyzer.



WARNING

If rack mounted, the Analyzer cannot support itself using a front rack mount kit alone. The instrument *must* be supported by a shelf or additional rails attached to the rack.



CAUTION

If the analyzer has been stored at less than 10 °C, allow the components to equalize to room temperature before starting the installation process.



CAUTION

Analyzer is for indoor use only and has an ingress protection rating of IPx-0. Analyzer is NOT protected against exposure to water including dripping, spraying, splashing or immersion.



CAUTION

When the analyzer is being integrated to an external system, the safety of that system is the responsibility of the assembler of that system.



CAUTION

During installation and operation, do not position the analyzer so that it is difficult to operate the disconnecting device.



NOTE

Take care to ensure that warm air is exhausted from an enclosure in which the analyzer is mounted.



CAUTION

It is imperative that the analyzer have adequate ventilation and/or cooling to maintain the ambient temperature below 35 °C when operating. Failure to provide adequate airflow, especially clearance at the front and rear panels, to ensure proper airflow and/or cooling to the analyzer will result in overheating of the analyzer causing a shutdown and potential damage. There should be 4" (10cm) of clearance in the front and back of the analyzer.

Thermal Specifications	Min	Max	Description
Ambient Operating Temperature	10 °C	35 °C	Worst-case environmental limits (unless otherwise specified)

7.2 Analyzer and Vacuum Setup

1. Remove the Analyzer and the External Vacuum Pump from their respective shipping containers.



This analyzer weighs 46 lbs. (20.9 kg). Use the technique outlined in the General Safety section on page 15 when lifting or moving the analyzer.

2. Place the Analyzer on a bench top or flat surface. Place the External Vacuum Pump near-by or on the floor. Don't push the analyzer into position yet, there are cables to be installed on the back panel.
3. Unpack the analyzer accessories. The Certificate of Compliance and USB drive should be stored in a safe place and may be required if you contact Picarro with questions about your analyzer.
4. Remove the caps from the analyzer's INLET and VACUUM connection ports. Save the caps; you should reinstall them when the analyzer is stored, moved or shipped.
5. Remove the cap from vacuum pump's inlet. Save the cap for later use. Reinstall the caps when the pump is stored, moved or shipped.
6. Connect one end of the vacuum hose to the pump: hand tighten the nut and then use an 11/16" wrench (not included) to make an additional turn of one flat (about 60 degrees).
7. Connect the analyzer to a power source using the supplied AC power cable.




Use the AC power cables supplied with the analyzer or a similarly rated cable. Check with Picarro technical support if you have questions about power cable replacement. An inadequately rated power cable can result in equipment damage.



Cords shall be RATED for the maximum current for the equipment and the cable used shall meet the requirements of IEC 60227 or IEC 60245. Cords certified or approved by a recognized testing authority are regarded as meeting this requirement. The connector type used should be: IEC320 C13.

8. If desired, attach a tube to the External Vacuum Pump exhaust port and direct to a safe place for venting the mixture of sample gases.
9. Connect the External Vacuum Pump to a power source using the other AC power cable.

Figure 5:
Analyzer Back
Panel

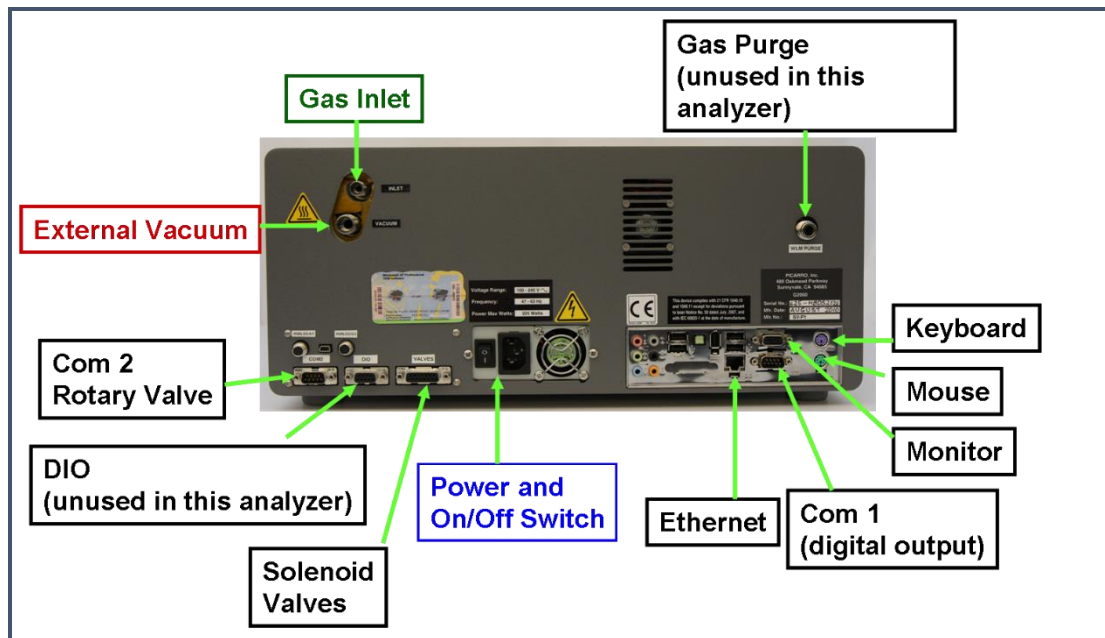

NOTE

The software to operate the instrument will start automatically after the operating system has loaded. The user interface will appear a few seconds after the instrument software starts. See “Startup Procedure” in Section 8.1.

7.3 Connecting to the Analyzer Inlet

Connect to the inlet of the analyzer using ¼" OD PTFE or PFA tubing using the supplied plastic ¼" PFA inlet nut and ferrules.



WARNING


When using compressed gases, follow all appropriate safety conventions, including use of eye protection, physical restraint of cylinders, etc.

Making a New Connection:

When using new tubing, follow these steps.

1. Place the two ferrules inside the nut as shown.

Figure 6:
Orientation of
ferrules and nut

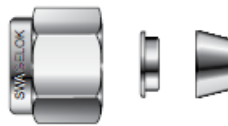

2. Loosely connect the nut to the INLET on the back panel of the analyzer about a full turn, being careful not to let the ferrules fall out.
3. Insert the tubing into the back of the nut, feeding it in as far as possible without deforming the tubing.
4. Hand tighten the nut.
5. Using a 5/8" wrench (not included), tighten the nut approximately seven flats (420 degrees).

Replacing a Connection

When reattaching tubing that already has a nut connected:

1. Inspect ferrules. If you see any damage, replace the ferrules and follow the directions above for making a new connection.
2. If there is no damage, hand tighten the connector to the analyzer's INLET.
3. Using a 5/8" wrench (not included), tighten the nut approximately one flat (60 degrees).

7.4 Setting Up a Monitor, Keyboard, and Mouse

A video monitor (not included), keyboard, and mouse are required for monitoring device operations, viewing or changing settings (including setting user permissions), or validating device performance.



NOTE

This device will operate under its default settings without any direct control of the internal software.

1. Connect a keyboard to one of the available USB ports.
2. Connect a mouse to one of the available USB ports.
3. Connect a monitor to one of the monitor ports. The analyzer will detect the connection and adjust the resolution to match the monitor.

4. Connect the monitor to a power source.
5. Turn on power to the monitor.

8. ANALYZER BASIC OPERATION

8.1 Startup Procedure

This section describes the steps to power on the system. These steps are sufficient for operating the instrument.



CAUTION

Always turn on the external pump before turning on the analyzer.

1. Switch on the external vacuum pump.
2. Switch on the main power to the analyzer.
3. Press the soft power button located on the front panel. The software will start automatically, and the analyzer will display the CRDS Data Viewer window. Data Viewer features are detailed in Section 9.
4. Wait approximately 45 minutes for the system to power, initialize, and stabilize.
5. Once the cavity pressure and operating temperature are at their proper setpoints and stable, the analyzer will start measuring.


As the instrument is starting up, it is normal for there to be a delay in reporting data. This can take several minutes depending on how long it takes for the internal temperature and pressure to reach their operating points, and it is normal during this time for some concentration readings to be negative or constant.



NOTE

The data selection pull-down menus will not be populated with the appropriate items until actual data is being reported in the graph. This is typically less than 30 minutes, but depending on ambient temperature, the analyzer can take up to 1 hour to stabilize.

Figure 7:
Picarro
Analyzer GUI
The heading shown may vary depending on the analyzer type.

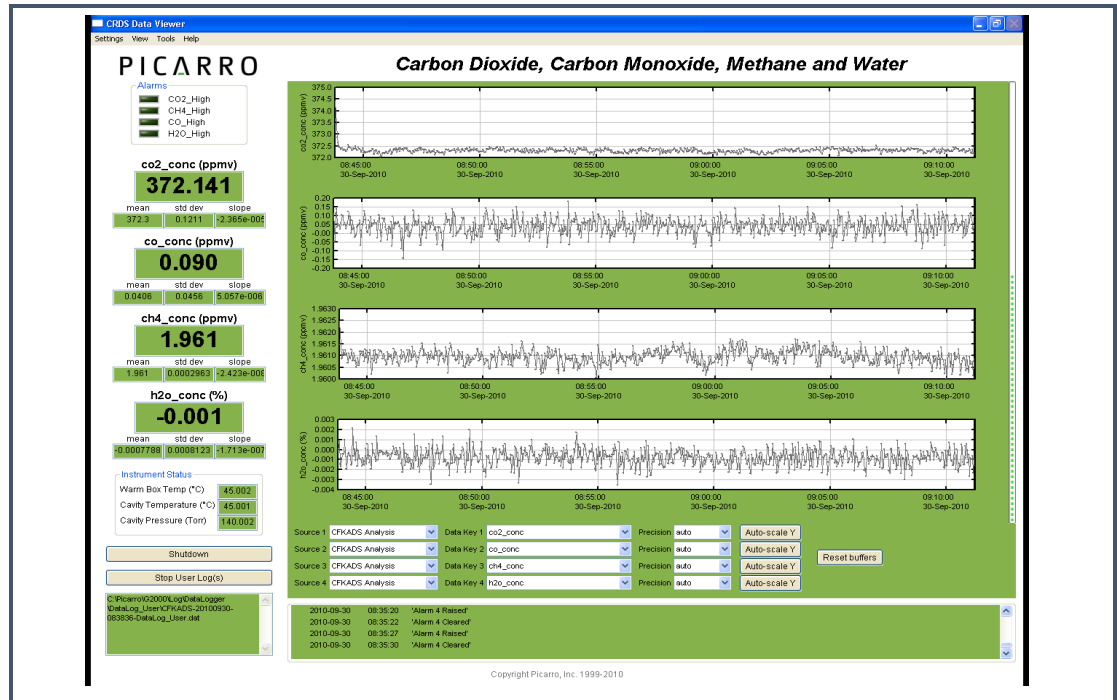


8.2 Desktop Icons and Folders

On the Windows™ desktop, there will be the following icons & folders:


- 1. Start Instrument:** When clicked, the analyzer will start measuring in the configuration that it was in when the software/analyzer was shut down.
- 2. Coordinator Launcher:** Depending on the system's configuration, the coordinator program may not be included. Clicking on this icon will lead you to a window that will allow you to choose the proper coordinator to operate the peripheral module that came with your analyzer.
- 3. Picarro Mode Switcher:** When clicked, you will be led to a window that will allow you to switch between various measurement modes. Most analyzer models are configured for one mode and may not include the Mode Switcher. If the analyzer has multiple modes, this allows the user to switch between them easily.

Figure 8:
Picarro Mode
Switcher

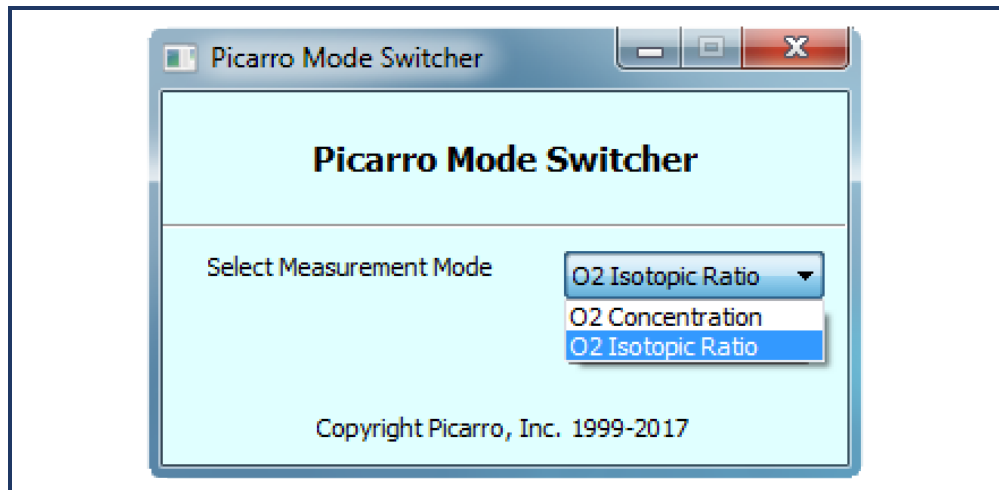

4. **Picarro Controller:** When clicked, you will be led to a useful diagnostic panel allowing the user to see the analyzer's internal temperatures, pressure, and spectroscopy in real time. This program has user-accessible functions, but it cannot change anything related to analyzer functionality. It is intended for diagnostic purposes only.
5. **Picarro Utilities Folders:**
 - a) **Data file Viewer:** When clicked, you will be led to a window that will allow you to convert between *.dat & H5 data files and to make various graphical representations of your data. The instructions on using the Data File Viewer software are described in the Appendix B
 - b) **Data Recal:** When clicked, you will be led to a window that will allow you to recalibrate your data based on known, certified data.
 - c) **Setup Tools:** When clicked, you will be led to a window that will allow you to edit various settings for your analyzer (See the Setup Tools section 10.3 of this manual for information about these settings).
6. **Diagnostic Folder**
 - a) **Stop Instrument:** When clicked, you will be led to a window that will allow you to turn off the analyzer in an emergency event. Upon clicking on this icon, the following window will pop up. Please see **Shutdown Procedure** on page 31 of this manual to shut down the analyzer in normal circumstances.

Figure 9:
Stop CRDS
Software

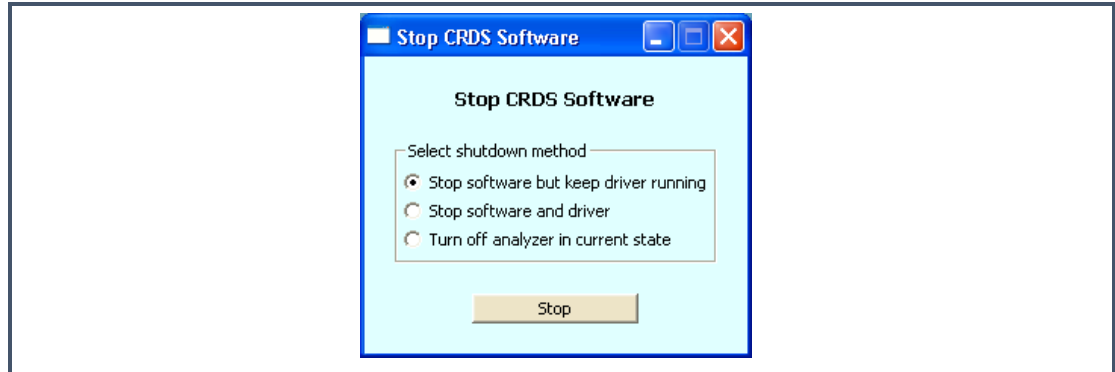

CAUTION

EXCEEDING GAS INLET PRESSURE OR TEMPERATURE SPECIFICATIONS COULD RESULT IN DAMAGE TO THE INSTRUMENT. In the case of higher input pressure or flow, configuring a sampling bypass manifold system is recommended.

Use a ‘tee’ at the gas inlet and either return the remainder to the main gas stream or exhaust appropriately.



CAUTION

Do not disconnect the AC power to the analyzer, vacuum line or the AC power to the External Vacuum Pump while analyzer is operating. Damage may be caused by current surges if power is applied while attaching or removing cables.

8.3 Shutdown Procedure

To shut down the analyzer using the GUI:

1. Click on the “Shutdown” button located on the left side of the Data Viewer window.
2. A window will pop-up prompting the user to confirm the shutdown. Once confirmed, the analyzer software and hardware will turn off.
3. Manually turn off the pump(s) and dry gas (only if your system requires it).



CAUTION

A flow of clean, relatively dry gas should always be directed to the instrument for several minutes prior to shutting down. Trapping a high-moisture content gas sample in the cavity can cause condensation damage to the mirrors as the instrument cools from its operating temperature.



CAUTION


Do not turn off the pump or disconnect the vacuum line while the instrument is operating.



CAUTION

FOR G2xxxx ANALYZERS: If you have trouble turning off the analyzer software, do not kill the process(es) in the task manager; rather, double-click on the "Stop Instrument" icon in the Diagnostics folder on your desktop.

Figure 10:
Shutdown
Window
The heading
shown may vary
depending on the
analyzer type

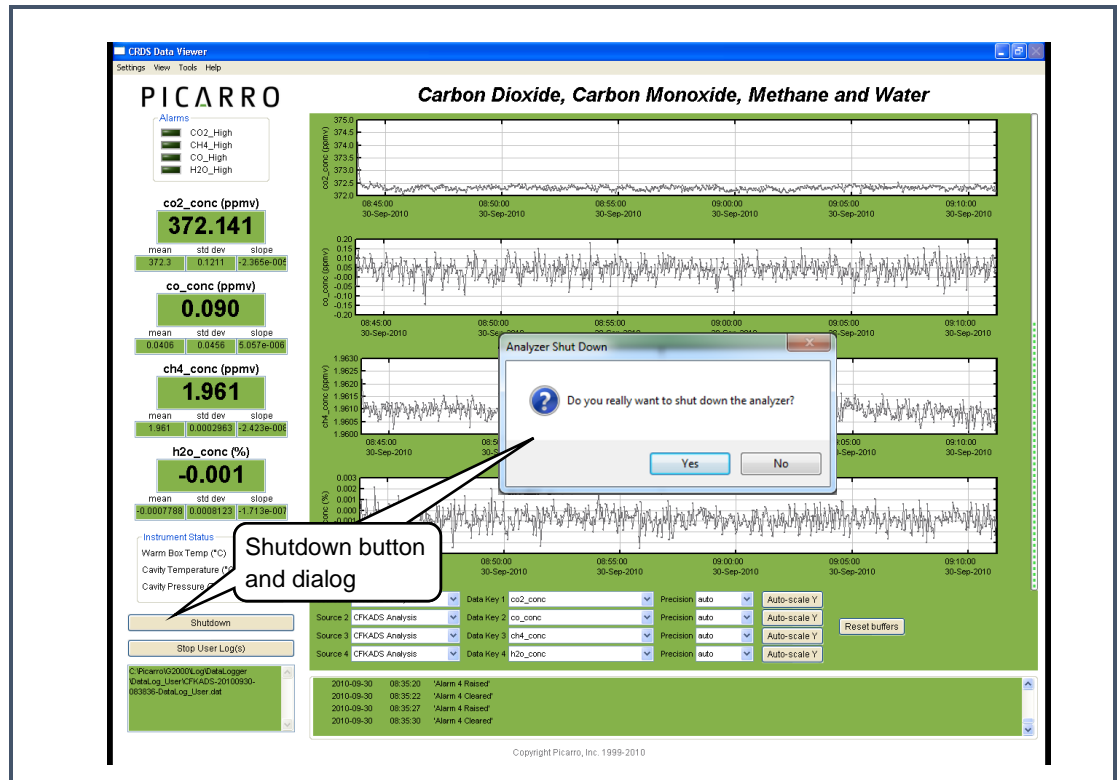

8.4 In case of an Electrical Power Outage

If the power to the analyzer is cut-off for any reason, the analyzer will cease operation. However, when the power is reapplied, the analyzer will restart automatically. The Picarro software tools will also properly close out previous files and open new files for data collection so that previously collected data, instrument diagnostics, and other parameters recorded up to the time of power outage are retained.


If short power outages will be a routine occurrence in your operating environment, Picarro recommends the use of an uninterrupted power supply to prevent damaging the hardware or software.

9. LIST OF GUI FUNCTIONS

9.1 GUI Overview

The features of the Picarro analyzer GUI are identified in Figure 11 and described in the following sections.

Figure 11:
Layout of
G2207-i
Analyzer
GUI

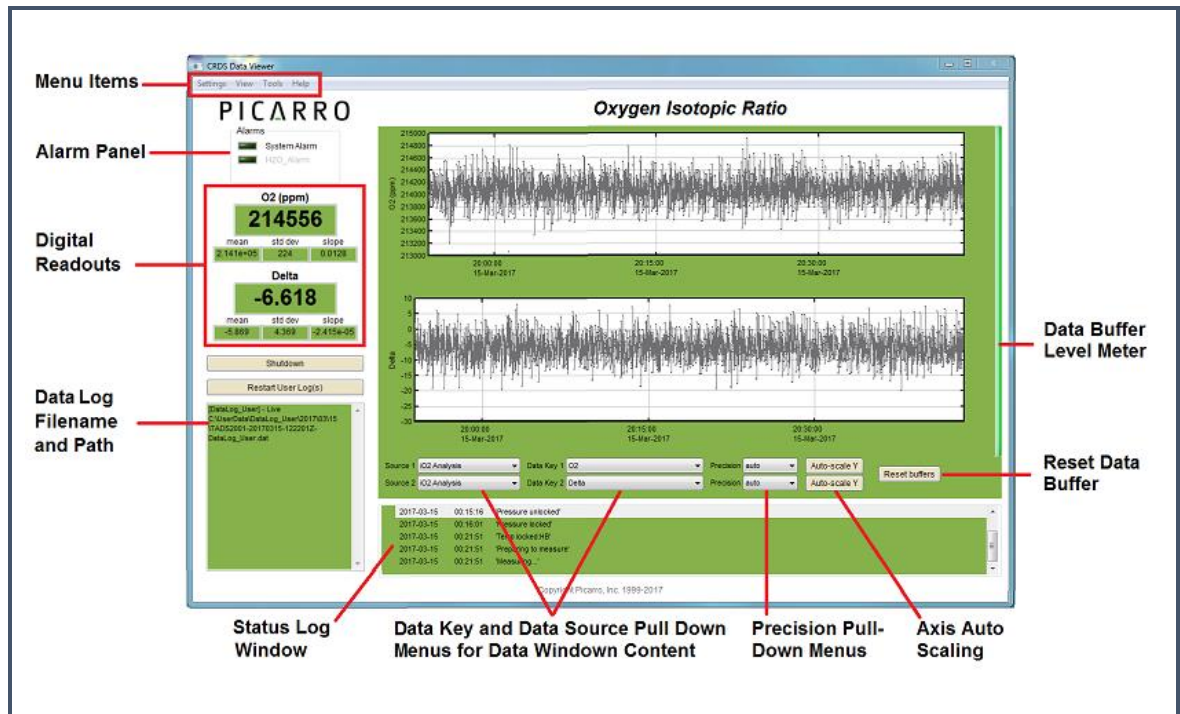

9.2 Settings, Tools, and Help Menus

Settings Menu

Left clicking on the Settings menu pulls down a menu that has one entry 'Change GUI Mode from Standard to Service.' This is the access point to a password protected service mode where additional operational and measurement parameters are displayed. Selecting and clicking on this entry opens the Cavity Ring-Down Spectrometer Controller. This is reserved for Picarro service operators only.

View Menu

This menu item has three entries:

1. Lock/Unlock time axis when zoomed: When locked, forces the two graphs to display the same time scale during zoom.

2. Show/hide statistics: Toggles the measurement statistics display, see **Digital Readouts** section on page **39** below.
3. Show/hide instrument status: Toggles the status of the instrument display.
4. Instrument Status below.

Tools Menu

This menu item has three entries:

1. **User Calibration:** Opens the user calibration window (default password is “picarro”). The password can be reset in the QuickGui.ini file in the instrument directory: “C:\Picarro\G2000\AppConfig\Config\QuickGUI” under the section: [Authorization] UserCalPassword = Picarro Show/Hide Valve
2. **Sequencer GUI:** Toggles the display of the external valve sequencer window.
3. The calibration slope and intercept can be entered and their effects immediately were seen in the data. Please refer to the Calibration Section on page 53.


Help Menu

1. **About:** Displays the version number of the instrument.

9.3 Alarm Panel

This panel is used to monitor the status of the internal instrument alarms. These indicators are gas concentration alarms, such as “NH3 Too High/Low” depending on instrument configuration. The gas concentration alarm LEDs are off (grayed) when the respective concentrations are below a certain value, and they are illuminated when the respective concentrations are above/below a certain value. To view the alarm set point, click on the LED and a dialog box will appear indicating the alarm setting and allow the user to enable it or change the setpoint:

Figure 12:
Alarm Setting
Window

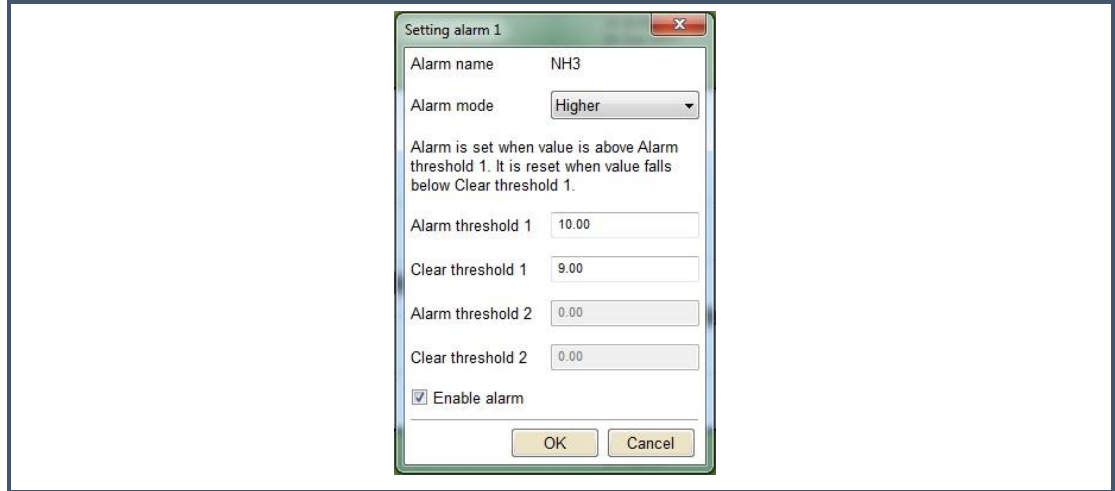


Type the value you wish to set the alarm to and press the “OK” button, or press “Cancel” if you do not wish to change the alarm value. If you do nothing, the dialog box will disappear and the alarm value will remain unchanged.

9.4 Instrument status

If these parameters are enabled through the ‘Show Instrument Status’ entry in the ‘View’ Menu on the main toolbar digital readouts for Warm Box temperature, Cavity Temperature and Cavity Pressure are displayed to the left of the main trend graphs.


Figure 13:
Instrument
Status

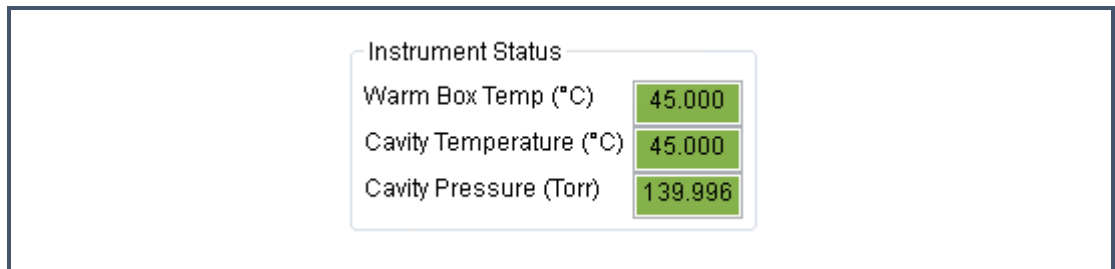

9.5 Data Log Filename and Path

The filename and path of the active data log are displayed in this pan. The indicator is grayed-out if there is no active data log (i.e., if a new data log has not been started using the *Start /Stop New Data Log Button*). A new file will be generated at midnight, which will be saved to the same location as the original log file.


Figure 14:
Data Log
Filename

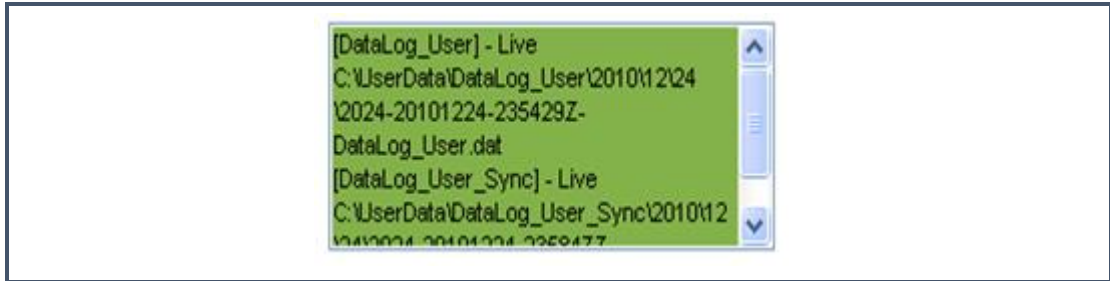

9.6 Data Window

The data window displays a graph of any stream of data vs. system time, with a format of hh: mm: ss. The user can select which data stream is displayed using combinations from the Data Source and Data Key pull-down menus. The precision displayed can be adjusted using the “Precision” menu and Auto-scaling of the ‘Y’ axis is also available.


Figure 15:
Data
Window

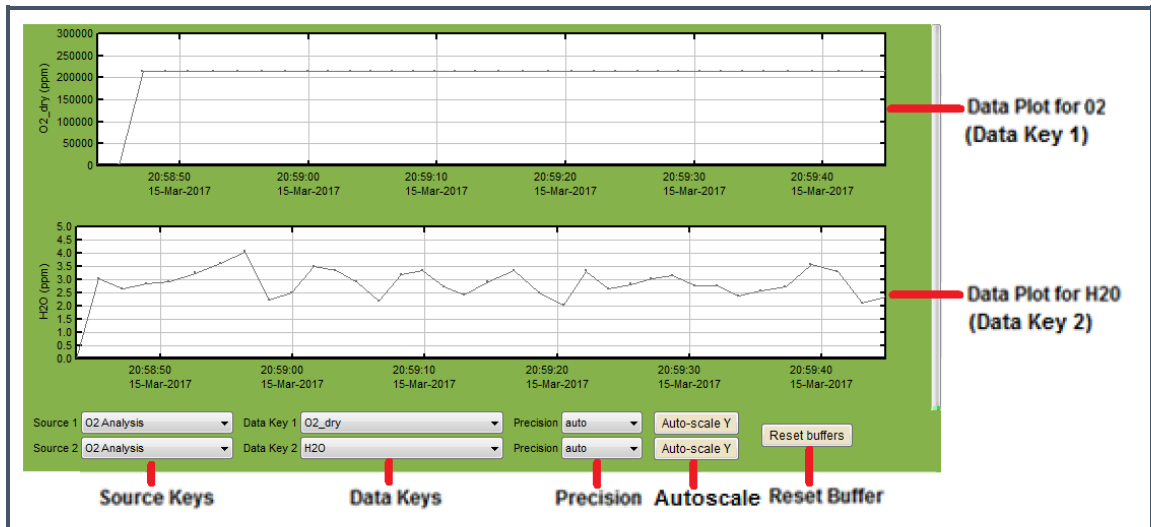

Data Source and Data Key Pull Down Menus

These two menus enable selection of the data stream that is viewed in the *data window*. Data streams available on the GUI are gas concentrations, if ‘*instrument Analysis*’ (where *instrument* represents the system installed) is selected, or if “sensors” is selected, the analyzer’s optical cavity pressure or temperature can be viewed as well as the nominal ambient temperature of the analyzer (“DAS temp”) and the temperature of the analyzer’s electronics chamber, indicated as “warm chamber temp.”


Figure 16:
Data
Source and
Data Key Pull-
Down Menus

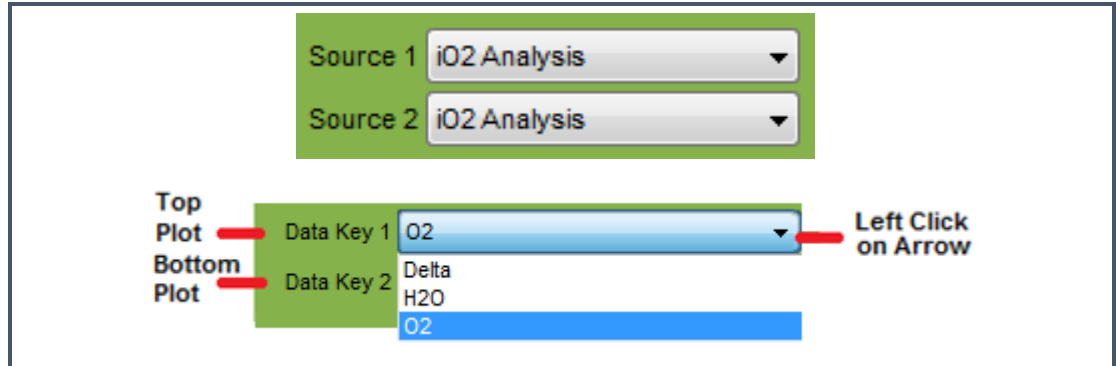

Precision Pulldown Menu

Click on this icon to select the precision displayed on the y-axis, between 0 and 4 digits of precision or “auto.” The currently selected precision is displayed during operation. This does not affect the precision of the saved data in the data log files or results files.


Figure 17:
Precision
Pulldown

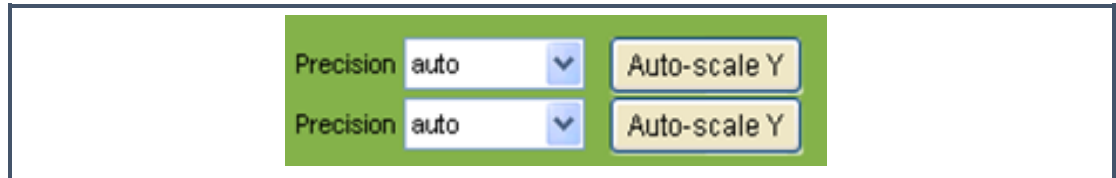

Data Buffer Level Meter

The meter to the right of the *Data Window* indicates how much of the internal memory of the GUI is used to retain historical data collected with the instrument. There is an internal limit of a finite number of points. Once that number of data points is collected, the buffer is full, and old data is removed from the buffer as new data is collected. This buffer affects only the data displayed in the data window, not the data stored in any files. This buffer is empty upon instrument startup, and can also be emptied by pressing the reset data buffer button in the lower-right-hand corner of the GUI.

Figure 18:
Data Buffer &
Reset Buffer

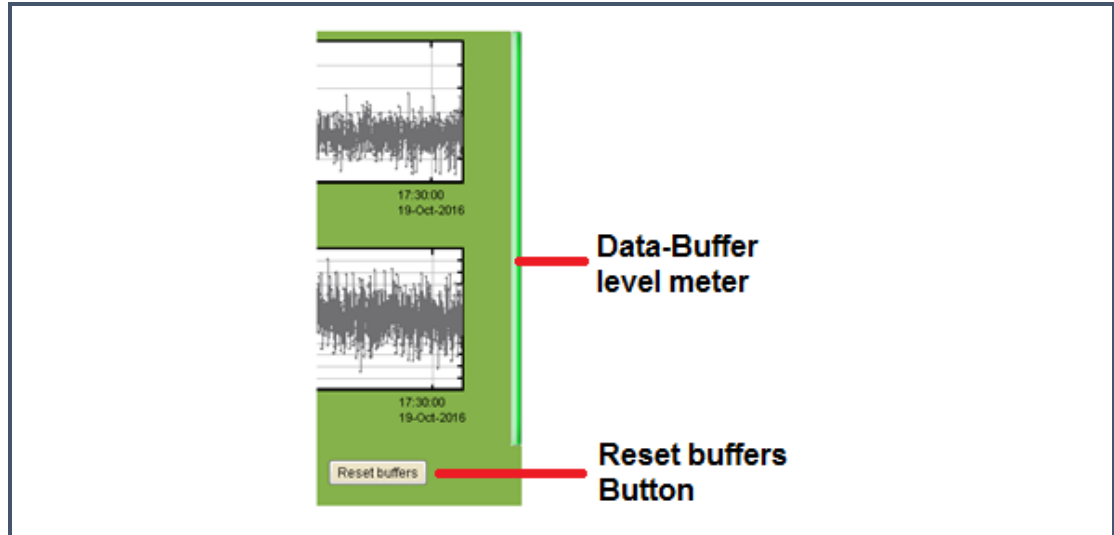


Reset Data Buffer Button

Press this button to clear the internal data buffer of the GUI (this clears the current data traces from the graphs). This has the effect of clearing all data in the data window. Pressing this button has no effect on any of the data log files stored by the instrument.

Graph Zooming

To zoom the graph, simply drag the magnifying glass over the section to be zoomed and click and hold the left mouse button. While holding down the left button, move the mouse to create a box that covers the region of interest. When the box is properly drawn, release the left button and boxed area will automatically scale to fill the data window. To zoom back out, double-click on the left button. To auto scale, the y-axis of either graph, use the auto-scale buttons below the graph. To lock or unlock the time axis of each graph during zooming, select that menu item in the 'View' menu.

Figure 19:
Graph
Zooming

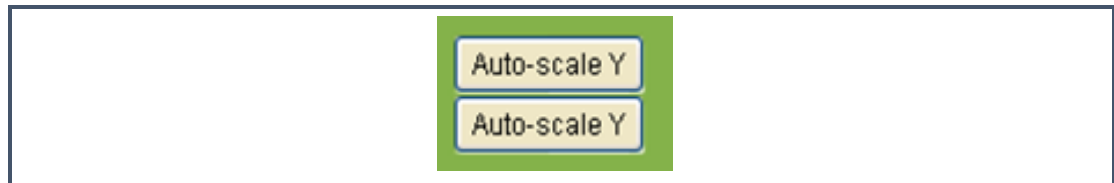


Digital Readouts

Displays the latest value recorded for the selected Data Key for each Data Window. Changing the Data Key modifies the Digital Readout as well as the Data Window view. If the 'Show Statistics' entry is enabled in the 'View' menu, the mean, standard deviation and slope of the data in the graph are

dynamically calculated and indicated below the digital concentration readout. These numbers change to reflect statistics of whatever data is in the data window.

Figure 20:
Digital
Readout

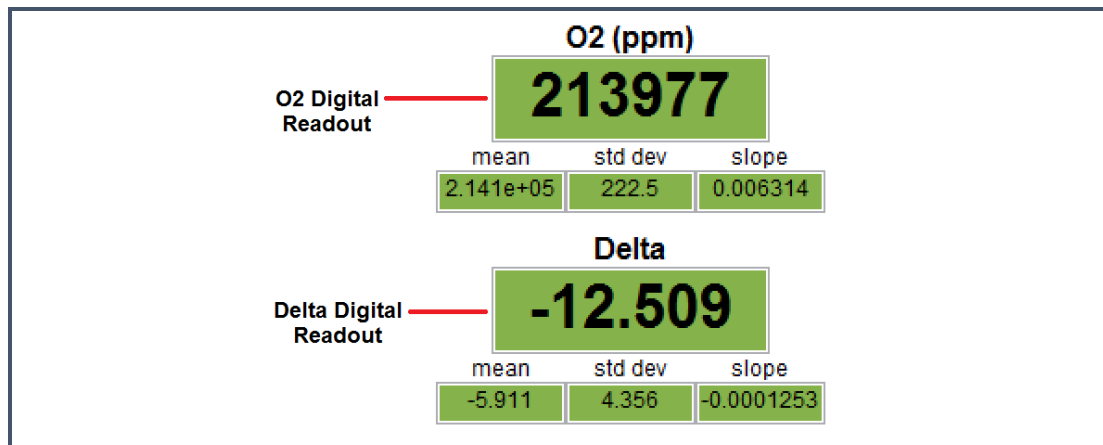

Start / Stop Data Log Button

The Analyzer automatically records all data collected on the instrument and saves it for later analysis. These files are called Data.dat files, which are described below in the section called **File Management**. In addition, the user can record a separate data log file. Press this button if you would like the instrument to start recording a separate data file. A dialog box will appear prompting you for a filename and location. Press this button again to stop recording the data file.

9.7 Status Log Window

This window displays instrument status messages, in the following form: “MM/DD/YYYY hh: mm: ss generic message text.” These messages include all messages sent to the DAS front panel display. The common status log messages are:

1. *Temperature Locked: WB (HB)*

The system waits for the warm box (“WB” – the temperature-controlled electronics and wavelength monitor chamber) to reach operating temperature. Similarly, the temperature of the hot box (“HB” – the temperature-controlled chamber containing the analyzer’s optical cavity and gas handling system) is stabilized. This is typically the longest step in the startup sequence. The duration of this step can range from 5 to 60 minutes, depending on the ambient temperature and how much time has elapsed since the last startup.

2. *Entering Measurement*

Spectral scanning has started. Concentration measurements will be available in approximately 30 seconds. The instrument will continue to scan and report concentration measurements until the instrument is shut down using the procedure below.

3. *Pressure Stabilizing/Locked*

The valve control system begins to allow flow through the analyzer and stabilizes the pressure inside the cavity.

4. *Measuring*

This is the normal mode of operation after startup has completed.

10. FILE MANAGEMENT

10.1 Overview

During operation, the Analyzer generates various ASCII-format text output files that are updated after each batch of concentration measurements is complete. For example, one of the user's output files is named:

CKFBDS##-YYYY-MM-DD-hmm- DataLog_User_Sync.dat

where "CKFBDS##" or similar is the instrument serial number.

The file name is generated from the instrument serial number, the date, and the time when the instrument was started. For example:

CKFBDS01-20070127-1029-DataLog_User.dat

- **CKFBDS01** = the instrument serial number
- **20070127** = the date, 1/27/2007, in format YYYY-MM-DD (to allow chronological sorting of data files).
- **1029** = the time the file was started, 10:29 am, formatted as hhmm using a 24-hour clock.

The raw user data is contained in folders in the directory:

C:\UserData\DataLog_User\[year][month][day]\hour.

This is data that has not been re-sampled to exact 0.1s time intervals.

There is a similar directory C:\UserData\DataLog_Sync\... which includes data that is evenly spaced in time at the data rate of the analyzer.

Data files are created every 15 minutes and stored for 90 days before they are automatically deleted. The file deletion frequency and details can be modified in the file:

C:\Picarro\G2000\AppConfig\Config\Archiver\Archiver.ini. This is explained below.

During data acquisition, the analyzer creates directories to store the data, based on the date the data were acquired. After each data file has been closed (every 15 minutes) it is moved to an archive directory, and a new file is started in the original location. The archive directory is C:\Picarro\G2000\Log\Archive\ and has subdirectories DataLog_Mailbox, DataLog_Private and DataLog_EventLogs with files arranged by year\month\day\hour.

There are complete data files which include additional information beyond the concentration data including parameters such as instrument temperatures and pressure, set points and spectroscopic information. This information is generally not useful to the user, but can be useful for diagnostic purposes and is stored in the directory
 C:\Picarro\G2000\Log\Archive\DataLog_Private
 \[year]\[month]\[day]\[hour].

For more information about how to include various columns of data from the DataLog_Private in the C:\UserData DataLog files, please contact Picarro.

The current data file can be found in C:\Picarro\G2000\Log\DataLogger\. Within this directory, there are subdirectories for DataLog_Private, DataLog_User, and DataLog_User_Sync. Similarly, the archive directory has subdirectories arranged by file type. The subdirectories are further organized by \[year]\[month]\[day]\[hour].

To keep the data files easy to manage and to limit the size of individual files and directories, the software automatically generates new files each time the instrument is powered up and at midnight each night. When new files are created at midnight, their file name will contain the new date and a time of 00:00.

For example, if the system was started at 10:29 am on 2/5/2019 it would create a file named:

20190205\CKFBDS01-20190205-1029-UserLog.dat.

Then at midnight, a new file will be created:

CKFBDS01-20190206-0000-UserLog.dat.

Specifies which directory to look in for files to delete.

`extension = data`

Specifies which files having what extension are to be deleted. If empty, it deletes all files.

`delete_time_hrs = 48`

Specifies how long to keep files prior to deletion.

10.3 Set Up Tool

In the desktop folder called Picarro Utilities, the Setup Tool can be launched by double-clicking on its icon. The tool allows the user to configure data file saving details, including which data elements are written to data files, digital data output (via serial port or TCP/IP), remote data delivery (via email), and general GUI properties.


Data Logger

The **Data Logger** tab allows the user to configure various data file saving details, including which data elements are written to data files.

- **Data Columns:** Controls which data elements are written to data files.
- **Hours of Each Log File:** Controls the size of each data document.
- **Enable Mailbox Archiving:** Enables archiving of data in the mailbox folder: *C:\Picarro\G2000\Log\Archive\DataLog_Mailbox*
- **Archived Directory Structure:** Specifies part of naming convention for data documents.
- **Total User Log Storage Size (GB):** Specifies the size of storage allowed for User Data (Recent Data).

After making the appropriate edits, click “Apply” to put changes into effect and then “Exit” to close the window.

Figure 21:
Data Logger
Setup Window

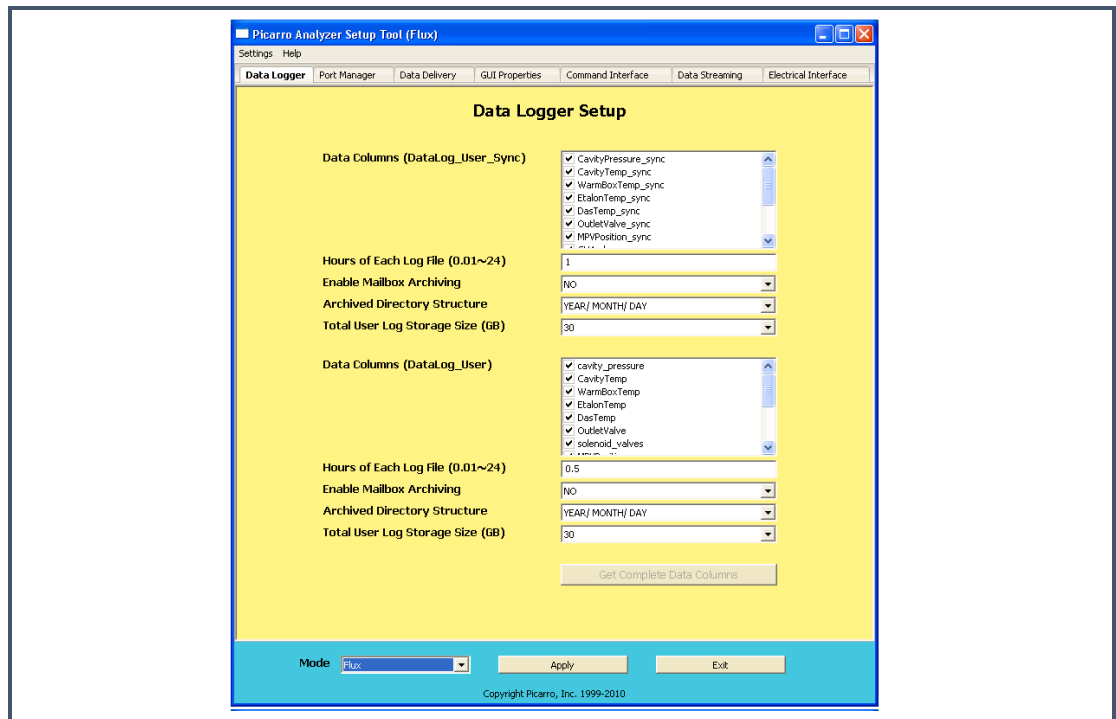


Serial/Socket Port Manager

The **Port Manager** tab allows you to control digital data output/Input via serial port or TCP/IP.

- **Data Streaming:** specify the port for data streaming (COM1/COM2/Off)
- **Valve Sequencer MPV:** specify the port for connecting to Valve Sequencer (COM1/COM2/Off)
- **Command Interface:** specify the port for Command Interface (COM1/COM2/ TCP/Off).

Make sure there are no COM port conflicts before clicking “Apply.” After making the appropriate edits, click “Apply” to put changes into effect and then “Exit” to close the window.

Figure 22:
Serial/Socket
Port Setup
Window

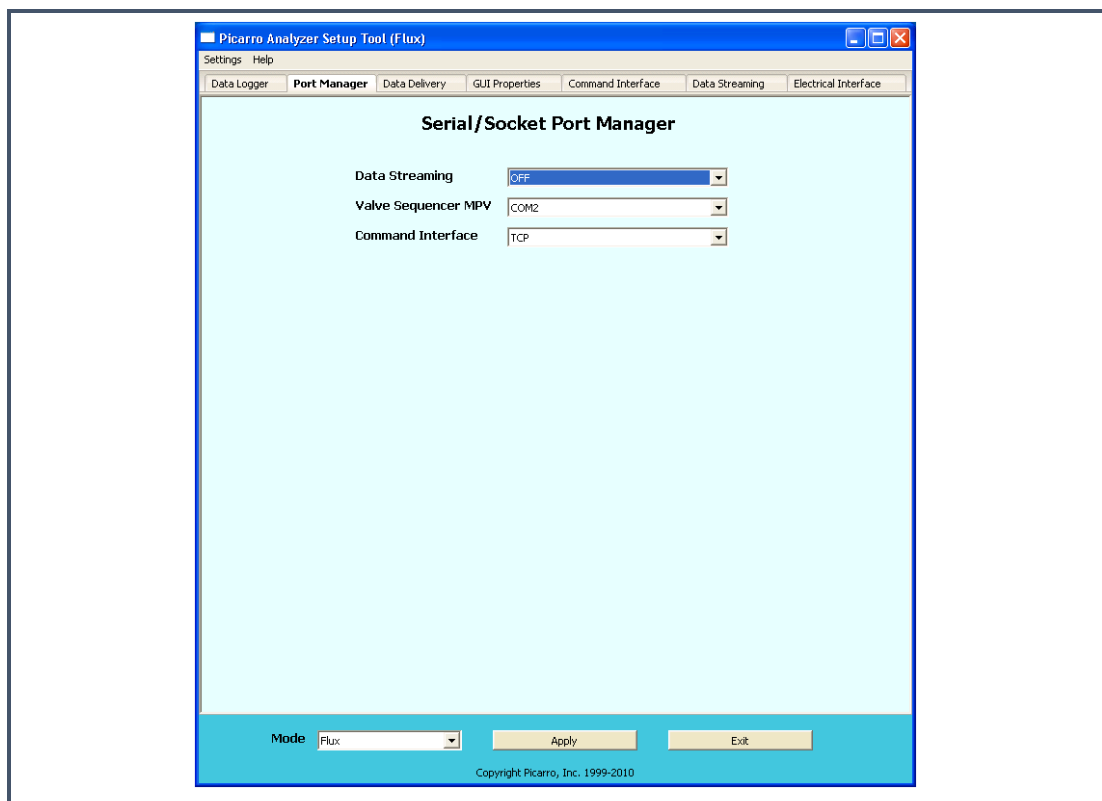


Data Delivery Setup

The **Data Delivery** tab allows the user to schedule remote data delivery via email.

- **Delivery Start Time:** Time of the day when data will be sent.
- **SSL:** Depending on the sender's email server, the sender can activate the Secure Sockets Layer (SSL).
- **Use Authentication:** Turning this on will require the receiver to provide a password and a username to access data. Set up the **password & Username** from this window.
- **From:** Sender's email
- **To:** Receiver's email.
- **Subject:** subject line of the email.
- **Data Directory:** Location of the data you want email.

After making the appropriate edits, click "Apply" to put changes into effect and then "Exit" to close the window.

Figure 23:
Data Delivery
Setup Window

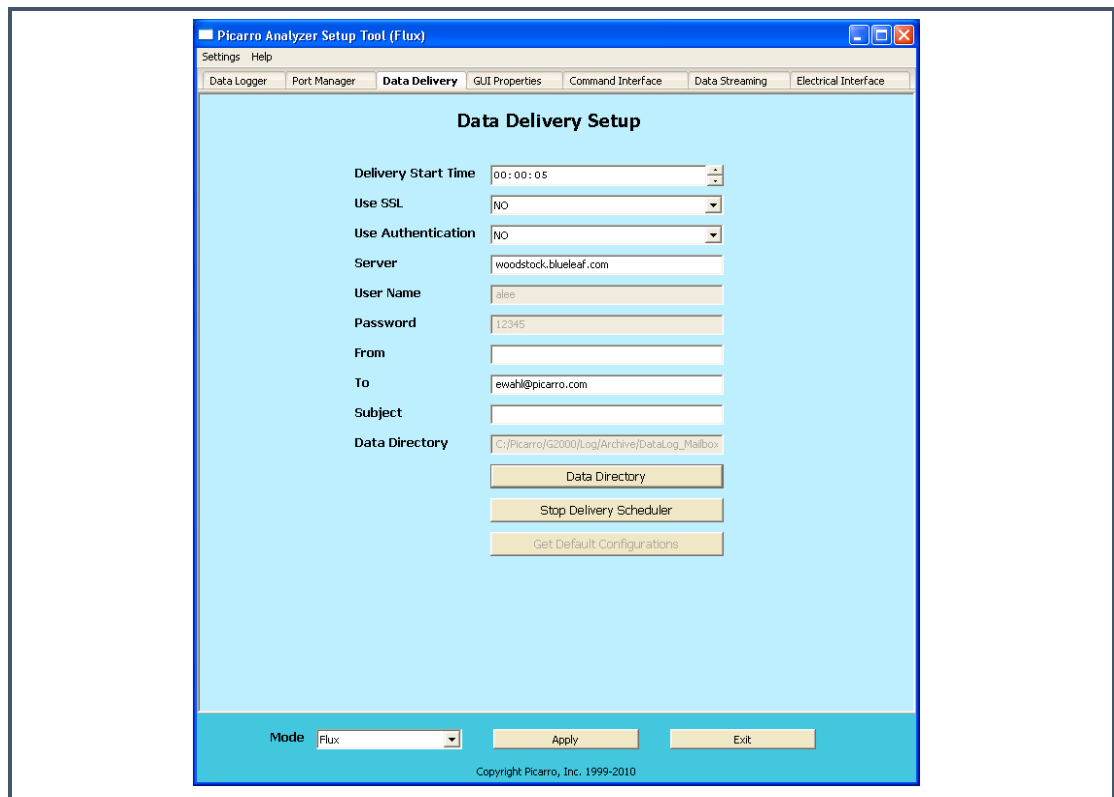

Editing Main GUI Properties


The **GUI Properties** tab allows you to set the number of line graphs visible on the main GUI.

In this tab, you can enable control of a connected Valve Sequencer from the main GUI:

1. Click on “Settings” of the “Setup Tool” window, and then “Switch to Service” mode.
2. Choose “Yes” next to “Enable Control Valve Sequencer” drop down menu on the “GUI Properties” tab.
3. Click “Apply” to put changes into effect and then “Exit” to close the window.

You should now be able to access the “Show/Hide Valve Sequencer GUI” menu from the main GUI under “Tools”

Figure 24:
GUI Properties
Window

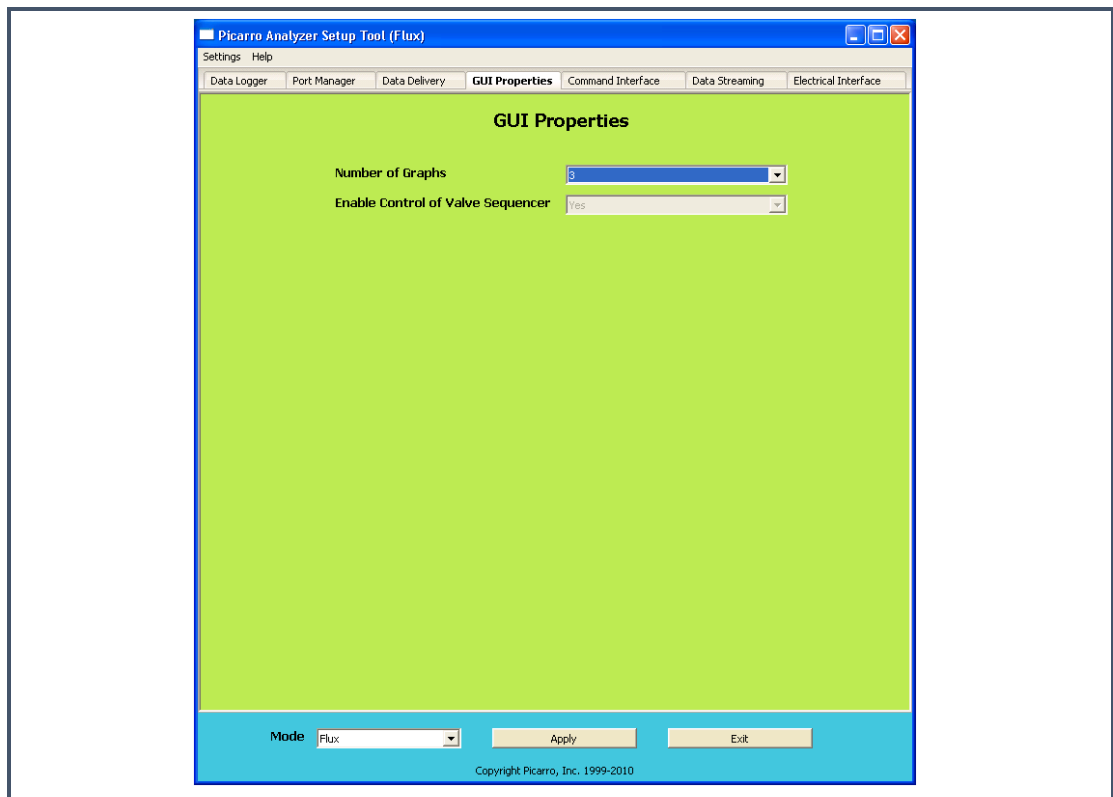


Specifying Digital Data Output for Command Interface

The **Command Interface** tab allows you to specify the data elements that are sent via COM port/TCP (specified in the Port Manager tab). Two types of data can be specified here: *DataLog_User* & *DataLog_User_Sync*.

After making the appropriate edits, click “Apply” to put changes into effect and then “Exit” to close the window.

Figure 25:
Command
Interface
Window

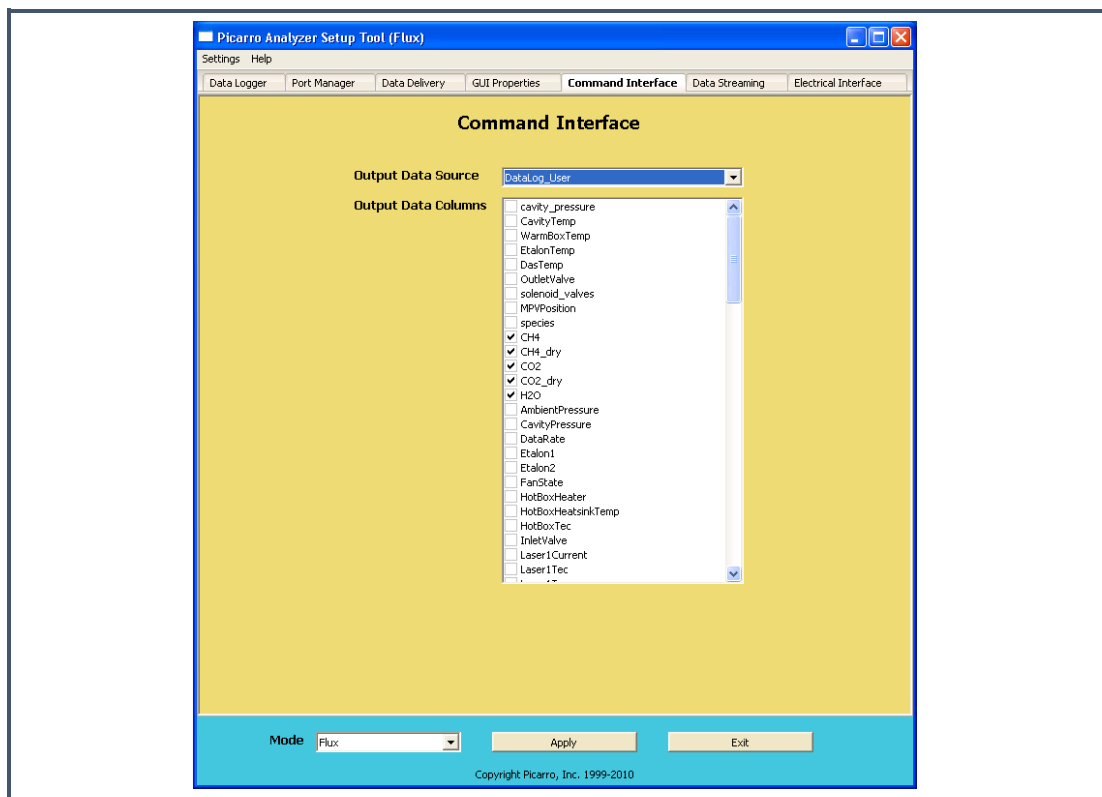


Specifying Digital Data Output for Data Streaming

The **Data Streaming** tab allows you to specify the data elements that you want to send via COM port (specified in the Port Manager tab). Two types of data can be specified here: *Datalog_User* & *DataLog_User_Sync*.

Data Streaming outputs data continuously, whereas the Command Interface needs commands to prompt data output.

After making the appropriate edits, click “Apply” to put changes into effect and then “Exit” to close the window.

Figure 26:
Data
Streaming
Window

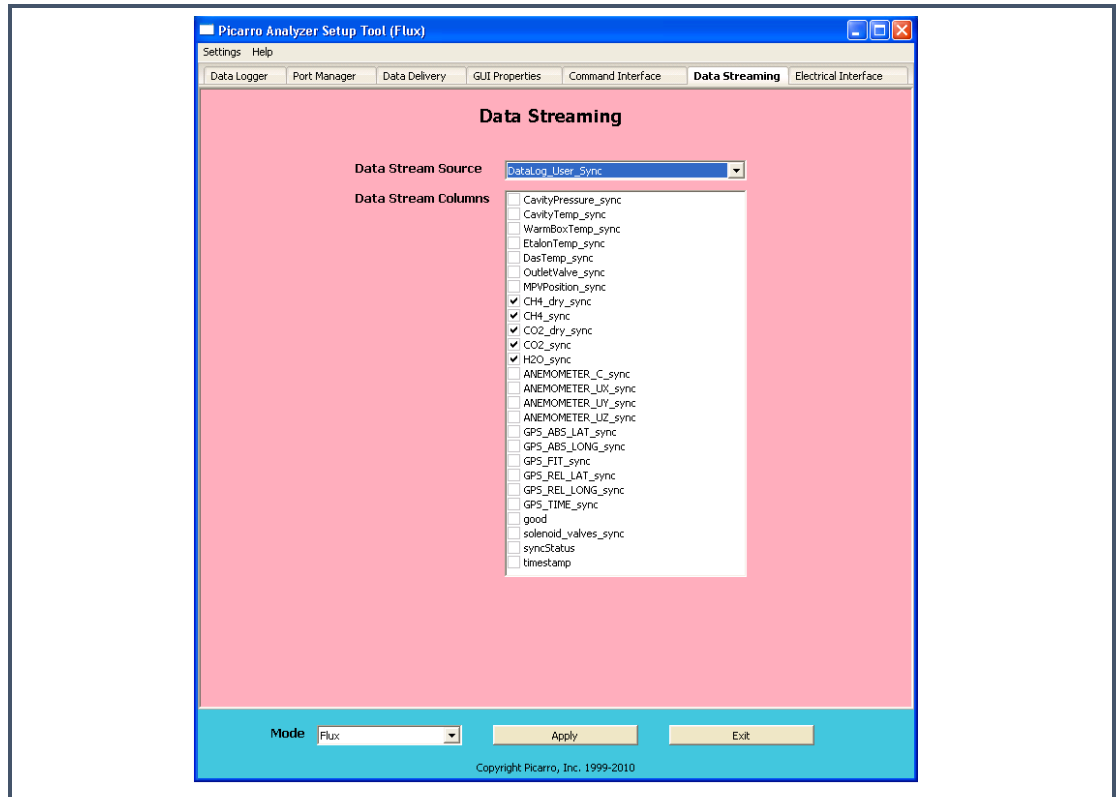

Customizing Analog Output Channels


The Picarro analyzer may be optionally configured with an Electrical Interface Card (EIC) that provides up to 8 analog signals available to the user for monitoring various measurements results and analyzer parameters.

Electrical Interface

The **Electrical Interface** tab allows you to customize each analog output channel. This tab will be disabled if your analyzer was not configured to work with an analog peripheral.

After making the appropriate edits, click “Apply” to put changes into effect and then “Exit” to close the window.

Figure 27:
Electrical
Interface
Window

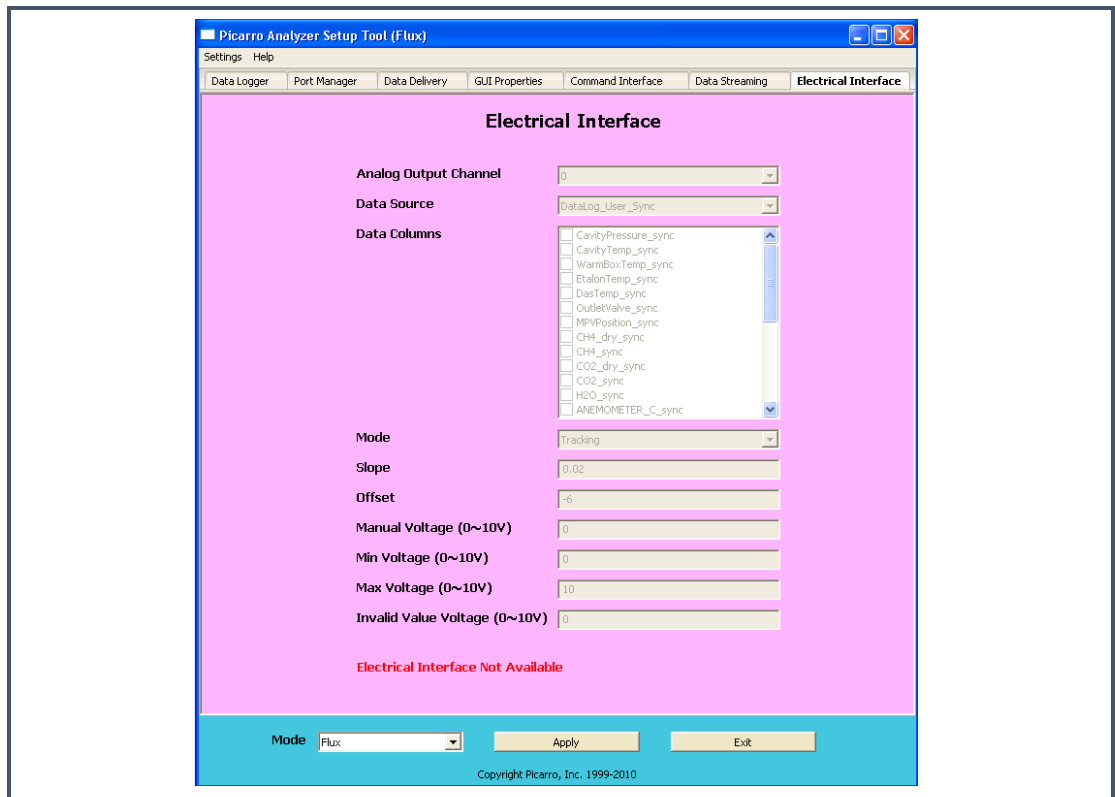

11. CALIBRATION

The two primary analytical measurements of the instrument are the concentration of NH_3 and H_2O . Periodic recalibration with standards of known concentration maintains the accuracy of these quantities. Using the Data Recal Software Utility (See instructions in Section 11.3: Calibration Software Utility Tools) enables the calibration constants to be tracked over time, thus enabling the user to follow system performance.

11.1 Slope and Offset

Every Picarro analyzer ships with a default user calibration value of 1 for slope and 0 for offset. While new user calibrations will likely result in changes to both, we recommend that only offset changes are applied during more frequent calibrations.

If a new slope calibration is to be applied, the user must be mindful of their experimental uncertainty. For example, a change in the slope calibration from 1 to 0.95 may not represent a change in the system's linearity, but may simply reflect the experimental uncertainty. If during your slope calibration, a new value of 0.9 or 0.85 is suggested by the Data Recal tool, please repeat the measurements of your standards. In most cases, such large changes in slope are a result of uncertainty and do not indicate a bad nonlinearity of the analyzer. If your analyzer continues to suggest a large change in slope, please consult with Picarro Support before applying a new calibration.



NOTE

If you are uncertain about the state of your current slope and offset values, they can be reset back to their user default values of 1 and 0, respectively. In your Picarro GUI, navigate to “Tools -> User Calibration -> (password: picarro) ->” Then change the slope and offset values of your parameter of interest to 1 and 0, respectively.

11.2 Calibration Methodology

To perform a calibration or verification of calibration, the user simply introduces the calibration standard into the analyzer for a period long enough for the analyzer to yield a stable measurement of that sample.

Calibration Setup

This section describes the connections from the analyzer to the gas tank.

- The pressure regulator at the outlet of the gas tank protects the analyzer from over-pressurizing. The pressure should be set to about 1 - 3 psi.
- The toggle valve allows rapid shutoff of the gas delivery.
- Tubing is connected to the male quick-disconnect connector provided with the analyzer.
- The male connector is inserted into the inlet port of the analyzer.

To replace the gas tank:

1. Turn off the main valve on the gas tank.
2. Disconnect the pressure regulator assembly from the tank.
3. Connect the pressure regulator assembly to the next gas tank.

Measurement Time for Each Standard

The measurement period for a calibration standard is dependent upon the recorded precision of the standard gas and performance characteristics of the analyzer, using the Allan standard deviation plot.

Measuring Multiple Gas Standards

When measuring multiple gas standards, the order of the gas standard is not important. However, make sure that you measure the gas standard slightly longer than the time you determined in the previous section because the usable/effective data set will be trimmed down.

Calibration Data Processing

1. Process the calibration results from the .dat file (see *File Management* on page 42) and calculate the average recorded value for each standard.
2. Plot these values versus the certified values from the gas supplier, and determine the linear relationship between the known calibration values and the analyzer's reported values. A linear best-fit equation can be calculated from the data.



NOTE

It is important to plot the analyzer’s reported concentration on the horizontal axis and the gas standards’ stated concentrations on the vertical axis.


The slope and intercept of the best-fit line through these points are the two values that are used to calibrate the analyzer. By determining the linear relationship between the known calibration values and the analyzer’s reported concentration values in this way, a calibration offset (slope and intercept) can be calculated to add a correction term to the analyzer’s factory or previous calibration.

Entering the Calibration Setting

Changing the analyzer’s calibration is intended to be done infrequently. Instead of recalibrating frequently to increase the accuracy of the data, users often just verify the calibration by measuring three or more gas standards and use the same regression procedure described here to calculate an offset by which to correct their data offline. Using the following equation in the graph below, this calculation would be accomplished point-by-point by calculating the corrected data “y” using the analyzer’s data “x” so that:

$$\text{Data corrected} = 0.9866 \cdot \text{Dataraw} + 5.268$$

Figure 28:
Linear
Calibration
Example

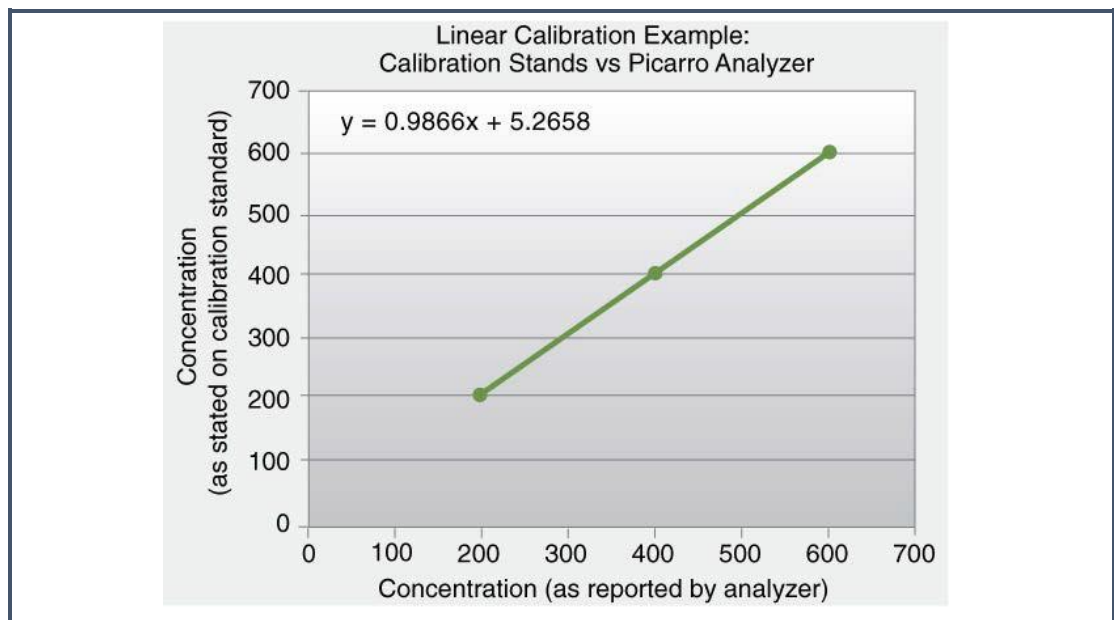

Calibration values are input into the analyzer by selecting **Tools > User Calibration**, and then entering the slope and intercept for each species.

**NOTE**

User Calibration is a password-protected form, and the default password is “picarro”. This password can be reset in the QuickGui.ini file as previously described.

The calibration will take effect immediately after clicking **OK**. To return to the factory calibration, simply set the slope to 1 and the intercept to 0 for each species.

11.3 Calibration Software Utility tools

Data Recal

The Data Recal software utility is designed to enable users to perform a routine concentrate calibration or calibration verification of Picarro’s gas concentration Analyzer via a user-friendly interface.

Data Recal allows the input of ten calibration points. Some of these points can be used for recalibration, while the remaining points can be used as quality control calibrants. In an ideal situation, three concentration certified standards should be analyzed on the system to generate CRDS-reported values for concentration. The standards with isotopic values spanning and encompassing the intended sample analysis should then be used to build a standard calibration curve. Parameters of this curve will then be used to correct the instrument readings so that they match standard values. The other standards will be simultaneously used for quality control, to verify other recalibrated delta values against their certified delta values.

Graphics User Interface (GUI)

The Picarro Data Recalibration (“Data Recal”) software can be found in the “Picarro Utilities” folder on the desktop. The Data Recal software utility can be launched by double-clicking on the Data Recal Icon in the folder. Shown below is a screenshot of the window that will open.

Figure 29:
Data Recal
Software Utility
GUI

The screenshot shows the 'Delta Recalibration' window with the following data and controls:

Used for Recal	Certified	CRDS Reported	Recalibrated
<input checked="" type="checkbox"/>	-35.6	-35.8	0.00000
<input checked="" type="checkbox"/>	8.6	7.6	0.00000
<input checked="" type="checkbox"/>	37.5	38.4	0.00000
<input type="checkbox"/>	13.6	13.2	0.00000
<input type="checkbox"/>	0.00000	0.00000	0.00000
<input type="checkbox"/>	0.00000	0.00000	0.00000
<input type="checkbox"/>	0.00000	0.00000	0.00000

	Current Calibration	New Calibration
Offset	1.87678	0.00000
Slope	0.54922	0.00000
R2		0.00000

Calibration Options: (dropdown menu)

Buttons:


Copyright Picarro, Inc. 1999-2009

The Delta Recal Software Utility comprises three sections:

1. A numerical input and selection section (shown below):

Certified and **CRDS reported** values are entered in the white boxes. Standards that will be used for instrument recalibration are then selected by checking the corresponding box in the first column, labeled “Used for Recal.”


Figure 30:
Recalibration
section of Data
Recal Software
Utility GUI

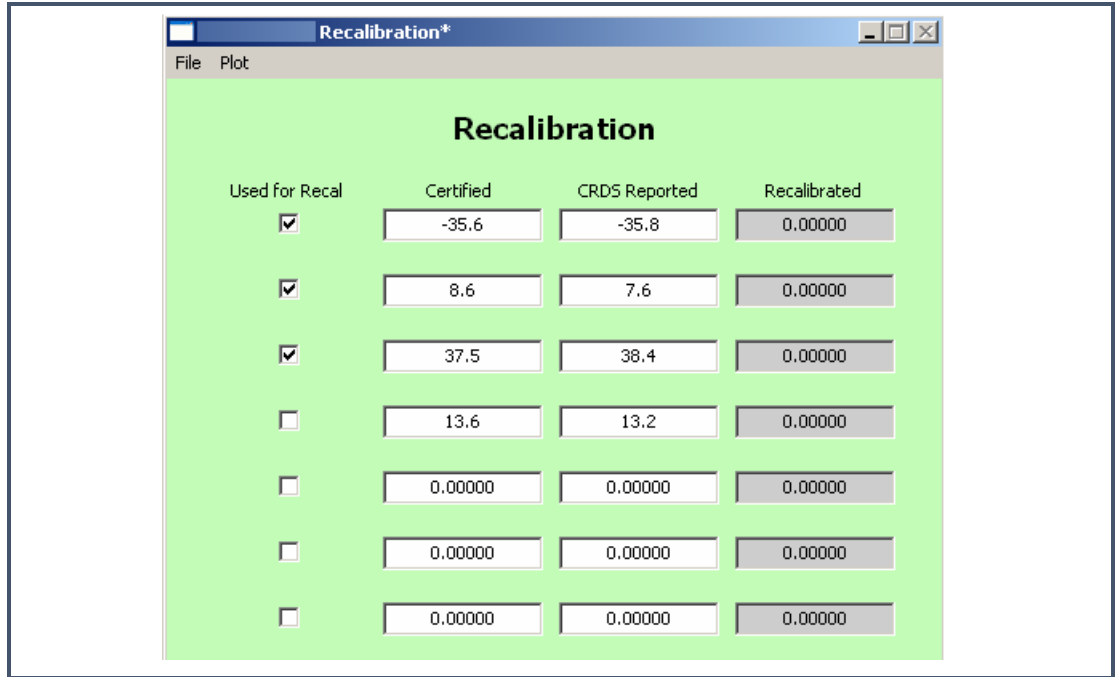

1. Calibration output section (shown below):

In the “Calibration Options” drop-down menu, you can select either an “offset” or an “offset + slope.” Once the new calibration parameters are calculated, the Data GUI will display the new values under the “New Calibration” column. The parameters that appear here will depend on your initial selection. When the offset+slope option is selected, the program also calculates a goodness-of-fit correlation coefficient (R^2).

Figure 31:
Calibration
output section
of Data Recal
Software Utility
GUI

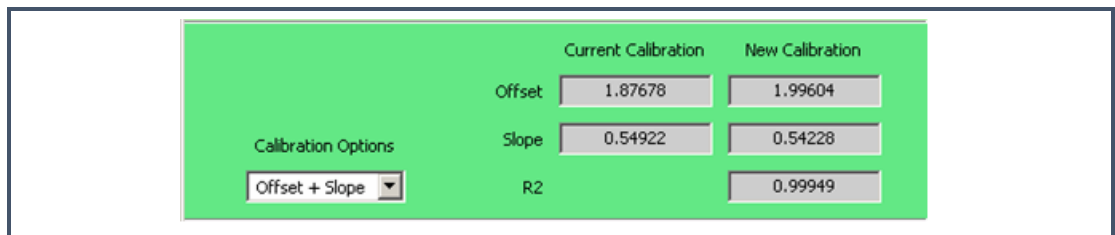


1. Action selection section (shown below):

In this section, the user can click on the compute button to calculate the new calibration parameter(s). The compute button will be grayed out until values have been entered in the **Certified** and **CRDS reported** columns, and at least one pair of **Certified** and **CRDS reported** values selected for the calibration. Once these values are entered or selected, the “Compute” button will be active. Click “Compute” to calculate the new calibration parameters. These

parameters will appear in the Calibration output section of the GUI, and an asterisk "*" will be displayed at the end of the window title line, indicating the new change.

Figure 32:
Action selection
section of the
Data Recal
Software Utility
GUI

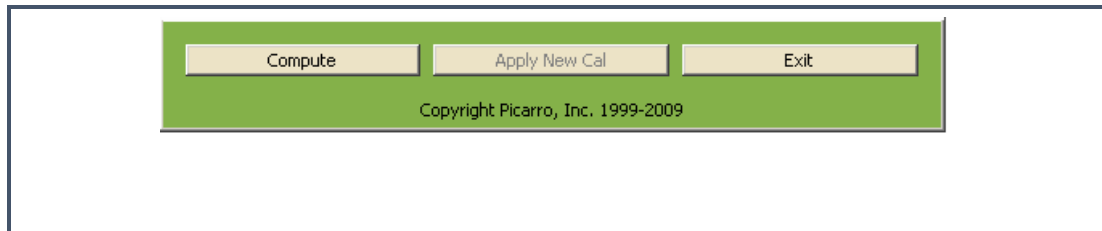

Once the calibration parameters are calculated, you can elect to graphically display the calibration and quality control (QC) standards in a graph. To do this, click "Plot" in the upper left corner of the window, and select "Plot Linear Fitting" from the drop-down menu, as shown below.

Figure 33:
Delta
Recalibration


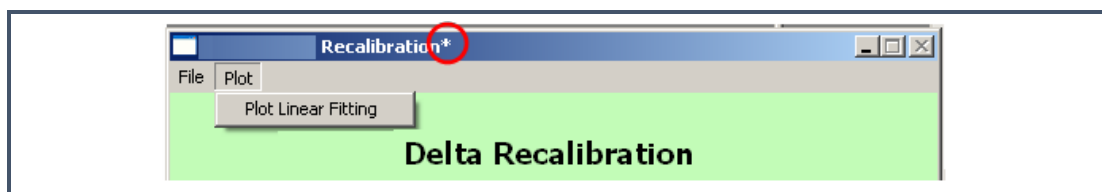


Figure 34:
Slope of Delta
Recalibration

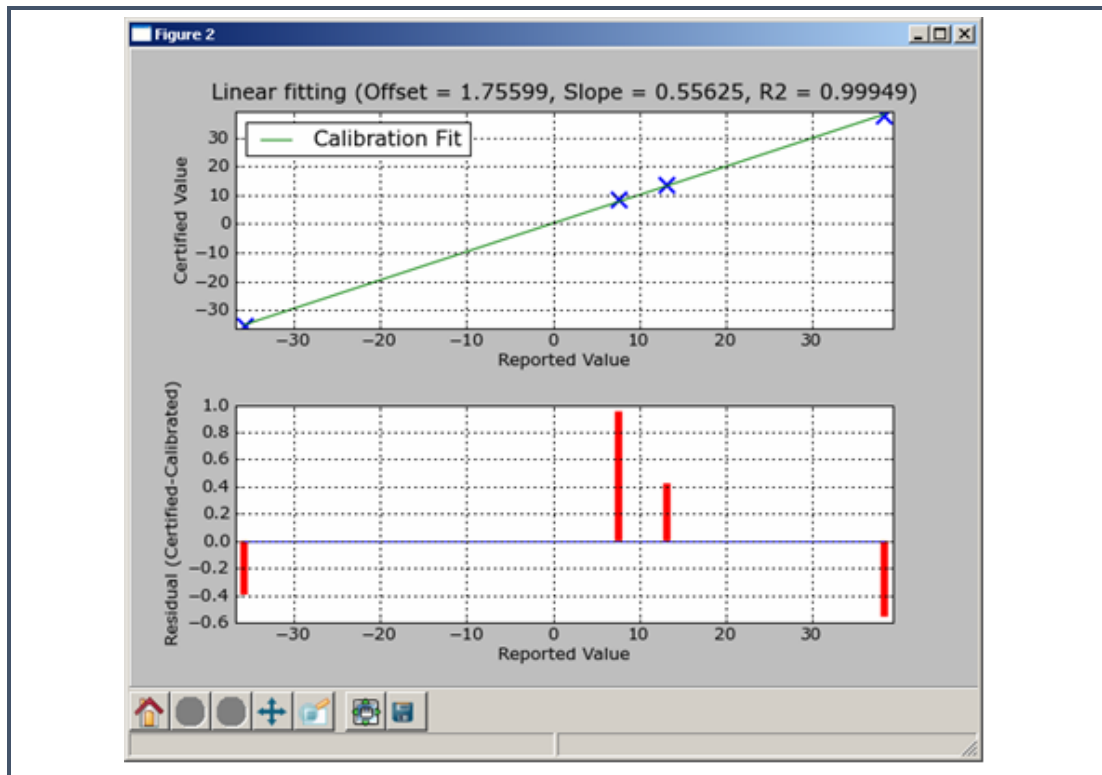

The plot that the utility will generate is shown below; this plot will display two graphs. The top graph presents measurements for the recalibrated standards (calibration and QC), as well as a line fit through those points. The bottom graph shows the residuals, which shows how far each measurement is from the best-fit line.


If you decide to accept the new calibration values based on this plot, click on the “Apply New Cal” button.

Figure 35:
Apply new
calibration slope
and intercept




Then you will be prompted to enter a Calibration Password. **The default Calibration Password is PICARRO.**

Figure 36:
User
Authorization
Dialog box

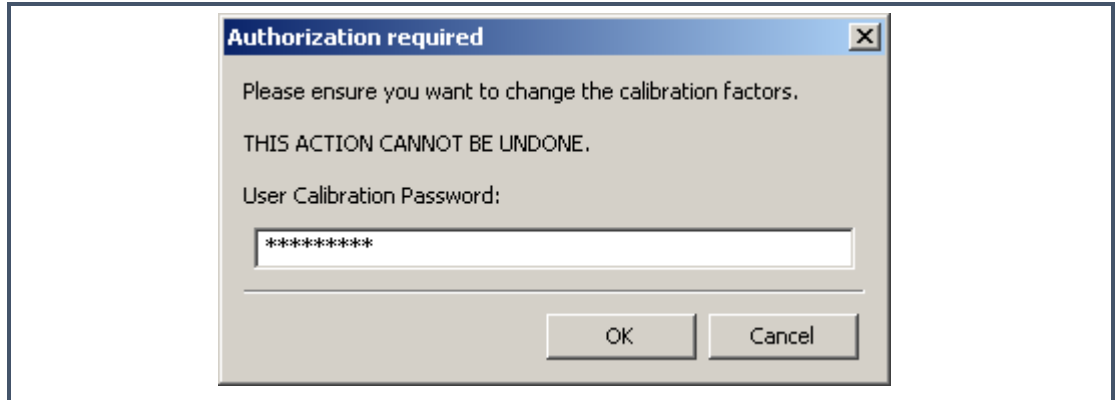

Next, a “Recalibration Confirmation” pop-up window will appear, displaying the new offset and slope (when applicable) values. This window will prompt you to confirm your choice to apply these values, to avoid any inadvertent mistakes.


Figure 37:
Calibration
Confirmation
Dialog box

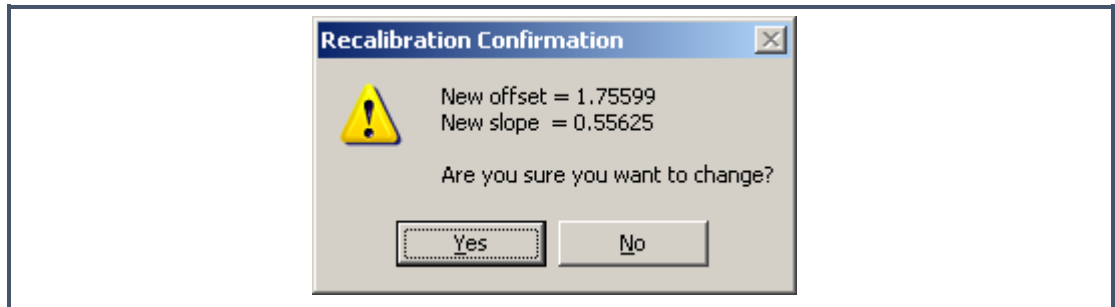

Please note that when you press the “Yes” button, the newly accepted calibration parameters will take effect immediately, without the need for the instrument main GUI to be restarted.


Once you accept the new calibration parameters, the Data Recal Software Utility automatically gives you the option to save the new recalibration file. Saving these files enables you to track the instrument recalibration history. Saved files can be reloaded by clicking “File” in the upper left corner of the window, and selecting “Load Recalibration File.”

Figure 38:
Delta
Recalibration
Load File

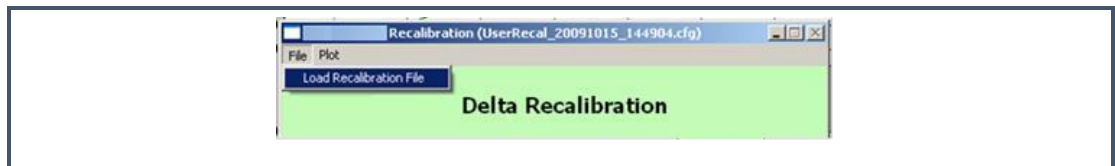

Next, select the file that you want to load.

Figure 39:
File Browser
Delta Recal File
Option

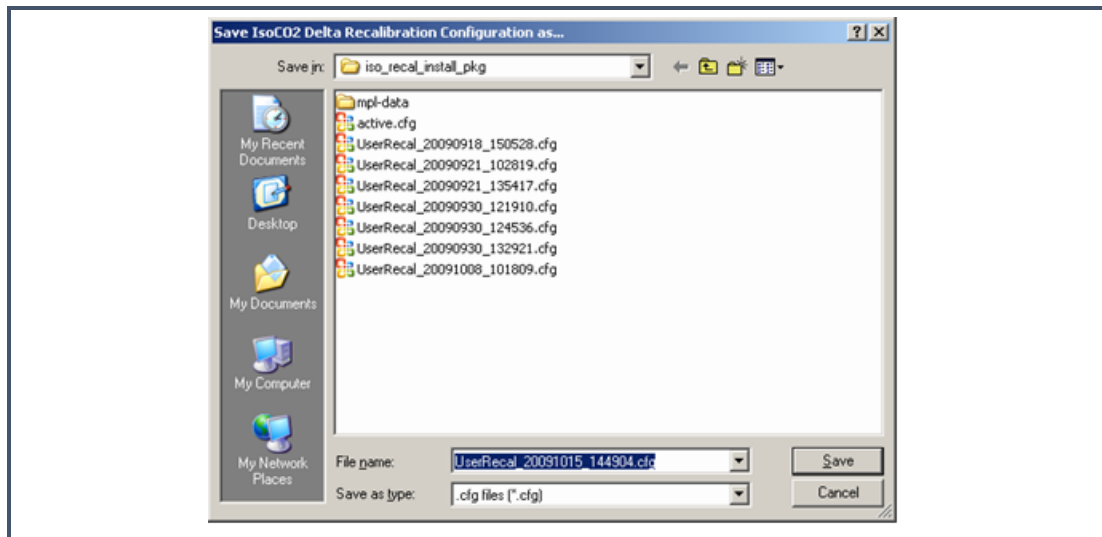

The recalibration file for this example will contain the following information:

Figure 40:
Delta Recal
Log File

-35.6	-35.8	-35.20437
8.6	7.6	7.64690
37.5	38.4	38.05748
cc	1.75599	0.55625
nc	1.87678	0.54922
		0.99949
option	Offset + Slope	


Rows 1 – 3: values from each standard selected for recalibration. The first column in these rows lists the certified values, the second column lists the measured values, and the third column lists the recalibrated values. The number of rows will depend on the number of standards selected to be used for recalibration.

Current Calibration (cc) Row: read from left to right, list the current offset and current slope value.

New Calibration (nc) Row: read from left to right, lists the new offset and new slope values. The third value in this row is the R^2 value, which is only displayed if the calibration option is “offset + slope.”


Options Row: The last row displays the calibration option selected. If you choose to exit the Recal GUI before accepting the new recalibration values, you will be prompted to confirm your choice through a pop-up window. This window contains a warning that continuing to exit will cause the new calibration data to be lost.

Figure 41:
Recalibration
Exit
Confirmation

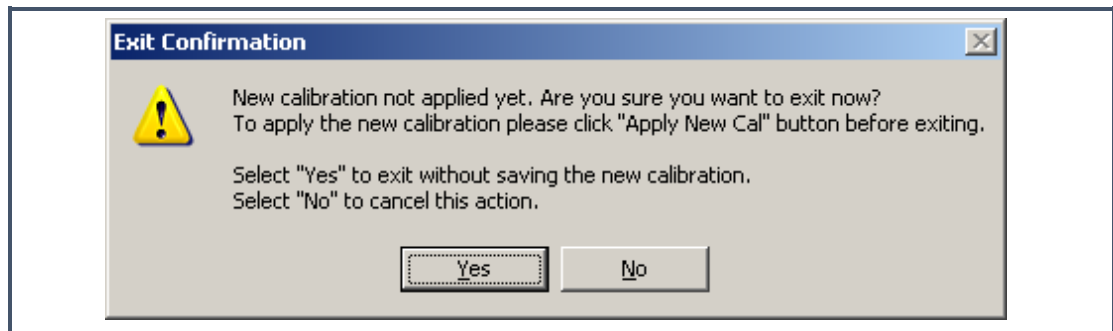

When you are finished with the calibration, you can exit the Delta Recal software utility by clicking on the “Exit” button at the bottom part of the Recal GUI.

Figure 42:
Calibration
Exit Button

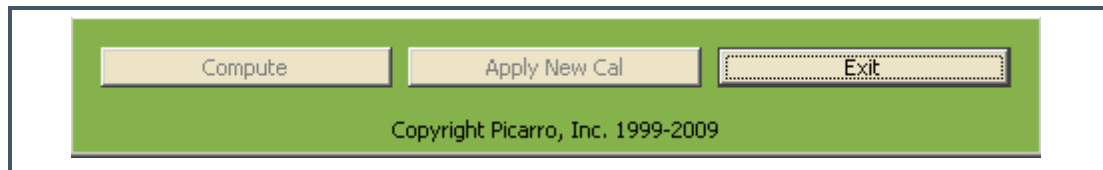

NOTE

The Recal GUI displays an error saying the new slope value is not acceptable when (1) at least two entries in the Certified or CRDS reported columns contain zero-value numbers, or (2) these entries are selected to be used for recalibration, or (3) the “offset + slope” calibration option is selected. This error occurs because the entered values will lead to an erroneous zero-slope value.

Figure 43:
Invalid New
Slope Entry



12. TROUBLESHOOTING

The following section lists problems that may be encountered during installation and operation of the analyzer. The corresponding step-by-step procedures provide resolution in most cases. If, after attempting these procedures, the problem remains unresolved, please contact Picarro Technical Support.

12.1 Power LED on Analyzer Does Not Illuminate

Context: Turning on the analyzer by momentarily depressing its front panel power switch should apply power. The green power LED is illuminated when it detects the correct power levels.

- a) Check that the AC power cord is attached and plugged into a working outlet.
- b) Check that the rear on-off switch near the AC power cord is in the on position.
- c) Press and hold the front panel power switch for at least 5 seconds as the analyzer may take several seconds to respond.

12.2 User Interface Program Does Not Start

Context: The computer may be configured to start the instrument and the associated user interface program automatically after it completes its boot-up sequence, or the program may be launched using the “Start instrument” icon on the desktop.

- a) Communications problems with the analyzer may occur if the analyzer fails to initialize correctly on power up. Should the analyzer initialization process not complete correctly, shut down the instrument by shutting down the Windows operating system on the control computer: use the Start menu, select the red Shutdown button and select “Shutdown” in the drop-down box under “What do you want the computer to do?”. Wait for the shutdown to complete normally and for the computer and analyzer to turn off completely. After a few seconds, restart the computer by momentarily depressing the power button.

Note: Do not simply restart Windows, since this does not cycle the power to the analyzer.

12.3 Sample Pressure Cannot Be Controlled to Appropriate Value for Concentration Measurements

Context: Under normal operation, the cavity pressure is automatically locked to the correct value by means of electronically controlled inlet and outlet valves. The message “Pressure Locked” on the front panel display and the user interface indicates that the cavity pressure is at the appropriate value. Should either of the messages “Pressure high” or “Pressure low” be displayed, the cavity pressure is out of its correct operating range.

- a) The “Pressure low” message indicates that there is insufficient gas available at the inlet of the analyzer. Check the inlet plumbing to the analyzer and ensure that the pressure at the inlet is within the specifications.
- b) The “Pressure high” message indicates that gas cannot be removed from the analyzer at a sufficient rate. Check the vacuum line between the analyzer and the power vacuum unit for leaks. Failure of the vacuum pump, injecting dilution gas at excessive pressure, or excessive pressure at the inlet can also cause this problem.

12.4 User Interface Program “Freezes”/Does Not Update Graphs as Data are Collected.

Context: The computer may become unresponsive causing the programs that control the analyzer to stop functioning. The computer and analyzer should be shut down and restarted.

- a) Resetting the computer and the instrument requires that the computer is shut down and restarted. If the computer responds to the mouse, a normal Windows shutdown may be carried out: use the Start menu, select the red Shutdown button and select “Shutdown” in the drop-down box under “What do you want the computer to do?” Wait for the shutdown to complete normally and for the computer and analyzer to turn off completely. After a few seconds, restart the computer by momentarily depressing the power button.
- b) If the computer does not respond to the mouse, hold down the power switch on the front panel for a few seconds until the computer and the instrument turn off. After another few seconds, restart the analyzer by momentarily depressing the power button.

13. SERVICE AND MAINTENANCE

The advanced, rugged design of the Picarro Analyzers provides stable, long-term operation with minimal service or maintenance. With the exception of the particulate filter, the analyzer is not user serviceable. Should it appear to malfunction, please refer to the Troubleshooting Guide or contact Picarro.

13.1 Particulate Filter Replacement

There are two in-line, sub-micron particulate filters before the measurement cavity. The first is user-replaceable and replacement filters can be purchased from Picarro and installed by the user.



CAUTION

Equipment Damage: It is important to **NEVER** remove the filter that is directly attached to the cavity. Only change the filter immediately following the inlet at the back of the analyzer. Refer to the filter replacement procedure in this document for further details.

The symptoms of a clogged filter can be analyzer reporting “pressure low” or there being no flow into the instrument, causing unusual measurements. Filters can become clogged after years of use in dirty environments. If liquid water is accidentally sucked into the inlet line, it will clog the filter and impede the flow (usually for a few days) until it evaporates. If this occurs, it is important to NOT turn off the analyzer or replace the filter until it is dry. The reason for this is that the increased humidity due to liquid water in the filter can cause condensation on the optics if the analyzer is allowed to cool from its operating temperature. Often, after the filter dries, the analyzer will begin functioning normally, and a filter replacement is not necessary.

Tools Required:


- 1.5mm hex driver
- 9/16” open-end wrench
- 5/8” open-end wrench
- 11/16” open-end wrench

Removing the Old Particulate Filter:

Move the analyzer to a clean work environment


2. Using the 1.5mm hex driver, remove analyzer top lid by removing six (3 per side) M3 x 6mm socket flathead screws.
3. Remove the piece of foam from around the input bulkhead by sliding it towards the back of the analyzer.

Figure 44:
Foam Removal

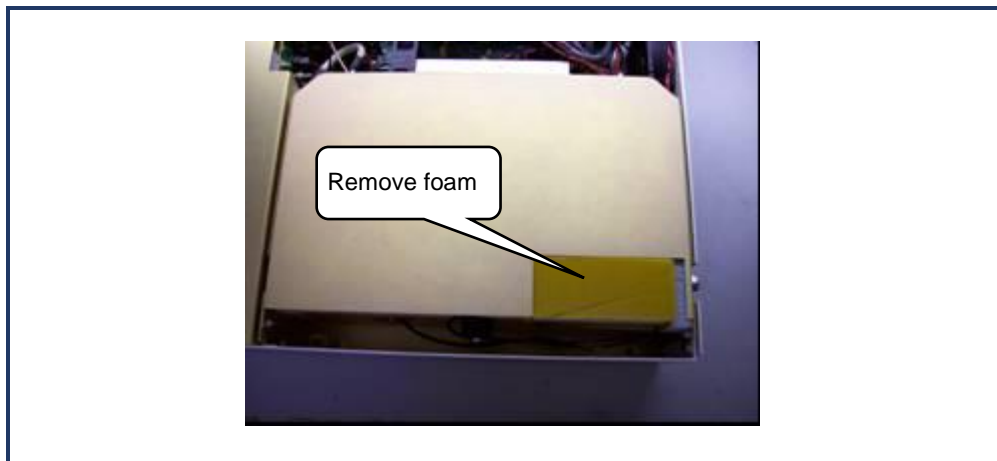


4. Using the 5/8" wrench, loosen the retaining nut on the input bulkhead (about 1 full turn should be enough).

Figure 45:
Loosen retaining nut.

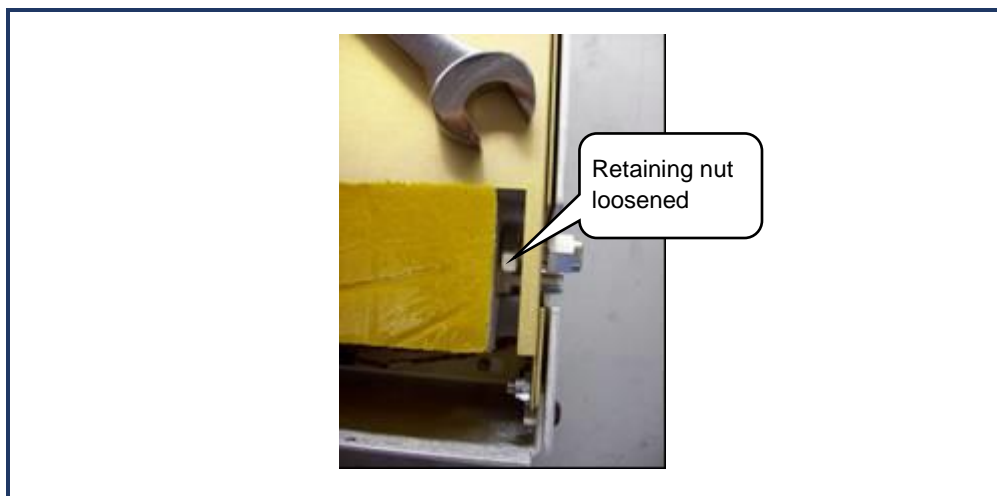


5. Slide the filter cover (with foam on top and side) towards the right side of the analyzer to remove it.
6. Using the 9/16" and 11/16" wrenches, disconnect the filter from the tube section near the front of the analyzer.

Figure 46:
Loosen nut at
filter output.

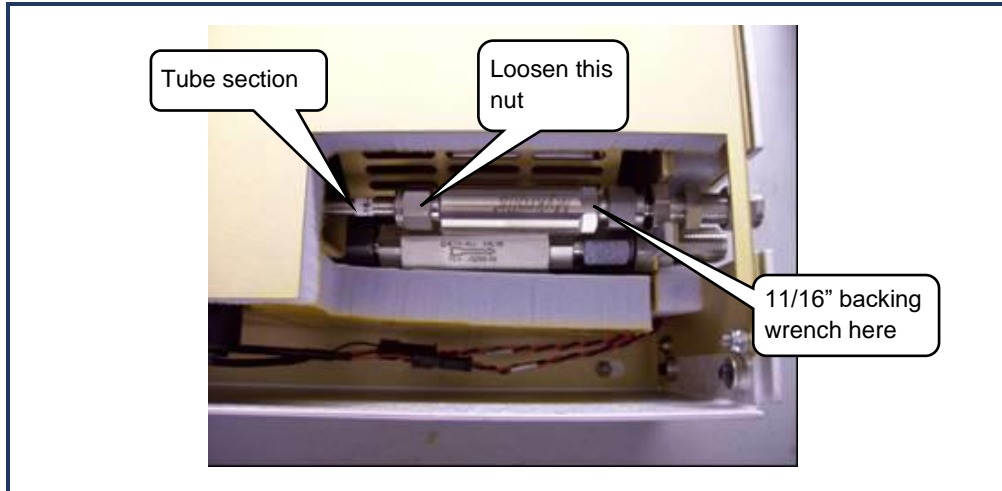


7. Slide the filter and bulkhead slightly towards the back of the analyzer and lift out.

Figure 47:
Filter and
bulkhead fitting
slid slightly
towards back
of DAS.

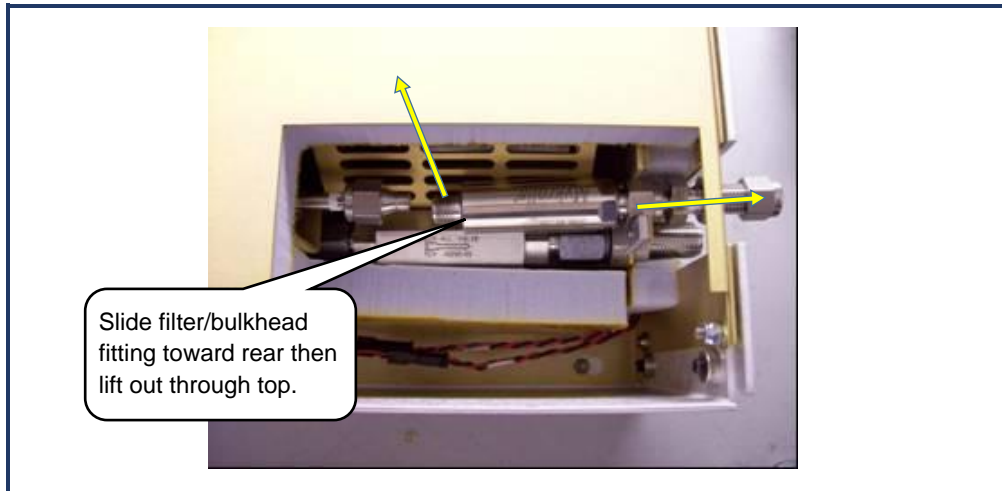

8. Using the 9/16" and 11/16" wrenches, disconnect the filter from bulkhead fitting.

Figure 48:
Filter removed
and separated
from bulkhead
fitting

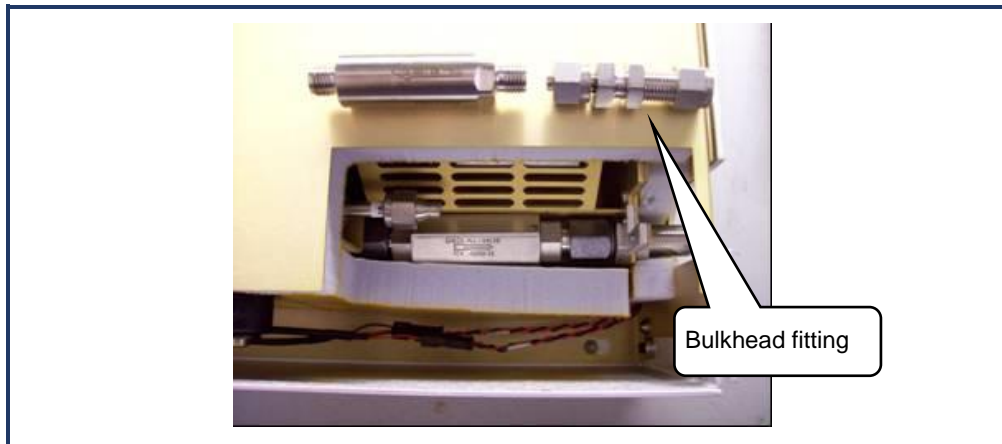

Installing the New Filter:



NOTE

When re-attaching 1/4" Swagelok fittings, the nut should be hand-tightened and then turned an additional 1/8 of a turn using a wrench.



CAUTION

Equipment Damage: To prevent cross-thread damage to the new filter during installation, carefully align the connector nuts with the filter before threading them on. Damaged filter threads can compromise the seal.




CAUTION

When installing the new filter, to ensure proper gas flow direction, make sure the arrow on the filter is pointing away from the back panel of the analyzer as shown in Figure 50.

1. Using the 9/16" and 11/16" wrenches, remove the filter from its packaging and attach it to the bulkhead fitting. The arrow on the filter needs to point away from the bulkhead fitting.
2. Using the 9/16" and 11/16" wrenches, reposition the filter and bulkhead fitting, and reattach to the tube section.

Figure 49:

New filter attached to bulkhead fitting. Note direction of gas flow.

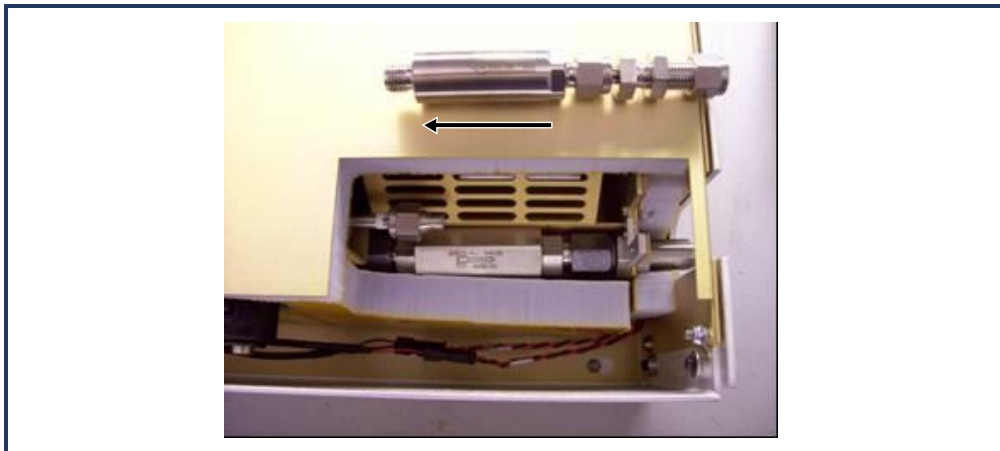

3. Using the 5/8" wrench, reposition the filter cover and tighten the retaining nut on the bulkhead fitting. The metal edge of the filter cover should be under the foam of the top of the enclosure.
4. Reposition the piece of foam around the input bulkhead fitting.
5. With the 1.5mm hex driver, reattach the analyzer top with 3 screws on each side.

13.2 Cleaning

Clean the outside of the analyzer with a clean dry cloth. Only certified service technicians should access or clean the inside of the analyzer.

14. TRANSPORTATION AND STORAGE

In the event that the instrument will be transported or stored, the following procedure can be used to prepare the instrument and repack it into the original carton.

Packing the Analyzer

1. Shutdown the instrument using the shutdown button, using the “prepare for shipment” option. Clean dry gas should be attached to the instrument prior to shutting down. This prevents condensation inside the system during storage or shipment.
2. Disconnect all tubing and electrical connections from the analyzer.
3. To prevent contamination and possible damage to the connector threads, place caps on all gas connections.
4. Place the analyzer in a plastic bag with a package of desiccant. Seal the bags with tape.
5. Pack the analyzer in the original shipping container ensuring that all the foam pieces are in place to protect the analyzer during shipping.




CAUTION

When shipping or relocating the analyzer, it is important to protect it from mechanical shocks. Failure to do so can compromise its performance. When shipping the analyzer, use its original packaging only.

15. LIMITED WARRANTY

Picarro, Inc. warrants its Products to be free from defects in material and workmanship and to perform in the manner and under the conditions specified in the Product specifications for twelve (12) months from shipment.

This warranty is the only warranty made by Picarro with respect to its Products and no person is authorized to bind Picarro for any obligations or liabilities beyond this warranty in connection with its Products. This warranty is made to the original Purchaser only, is non-transferable and may only be modified or amended by a written instrument signed by a duly authorized officer of Picarro. Sub-systems manufactured by other firms, but integrated into Picarro Products, are covered by the original manufacturer's warranty and Picarro makes no warranty, express or implied, regarding such sub-systems. Products or parts thereof which are replaced or repaired under this warranty are warranted only for the remaining, un-expired portion of the original warranty period applicable to the specific Product replaced or repaired.

	<p>DISCLAIMER</p> <p>THE FOREGOING WARRANTY IS EXCLUSIVE AND IN LIEU OF ALL OTHER WARRANTIES WHETHER WRITTEN, ORAL OR IMPLIED, AND SHALL BE THE PURCHASER'S SOLE REMEDY AND PICARRO'S SOLE LIABILITY IN CONTRACT OR OTHERWISE FOR THE PRODUCT. PICARRO EXPRESSLY DISCLAIMS ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE.</p>
---	---

The Purchaser's exclusive remedy with respect to any defective Product shall be to have Picarro repair or replace such defective Product or credit the Purchaser's account, whichever Picarro may elect in its sole discretion. If it is found that any Product has been returned which is not defective, the Purchaser will be notified and such Product returned at the Purchaser's expense. In addition, a charge for testing and examination may, at Picarro's sole discretion, be made on any Product so returned.

These remedies are available only if: i) Picarro is notified in writing by the Purchaser promptly upon discovery of a Product defect, and in any event within the warranty period; ii) Picarro's examination of such Product discloses to Picarro's satisfaction that such defects actually exist and the Product has not been repaired, worked on, altered by persons not authorized by Picarro, subject to misuse, negligence or accident, or connected, installed, used or adjusted otherwise than in accordance with the instructions furnished by Picarro.

The following warranty conditions shall apply to all Picarro, Inc. products unless amended by a written instrument signed by a duly authorized officer of Picarro:

ADJUSTMENT – No electrical, mechanical or optical adjustments to the product(s) are permitted.

PARTS AND LABOR - New or factory-built replacements for defective parts will be supplied for twelve (12) months from date of shipment of the product. Replacement parts are warranted for the remaining portion of the original warranty period. There will be no charge for repair of products under warranty where the repair work is done by Picarro, Inc.

NOT COVERED BY THE WARRANTY – Damage to any optical surface from improper handling or cleaning procedures. This applies specifically to those items subjected to excess laser radiation, contaminated environments, extreme temperature or abrasive cleaning. Damage due to ESD, abuse, misuse, improper installation or application, alteration, accident, negligence in use, improper storage, transportation or handling. No warranty shall apply where the original equipment identifications have been removed, defaced, altered or where there is any evidence of alterations, adjustments, removal of protective outer enclosure, any attempt to repair the product by unauthorized personnel or with parts other than those provided by Picarro, Inc.

DAMAGE IN SHIPMENT - Your analyzer should be inspected and tested as soon as it is received. The product is packaged for safe delivery. If the product is damaged in any way, you should immediately file a claim with the carrier or, if insured separately, with the insurance company. Picarro, Inc. will not be responsible for damage sustained in shipment. All Picarro products are F.O.B. origin, shipped from the Picarro factory or Picarro distributor. The price of all Products, unless otherwise specifically stated, is Ex- Works, Sunnyvale, CA as defined by Incoterms, 2001. The cost of normal packaging for shipment is included in the invoiced price. Where Buyer specifies special packaging, a charge will be made to cover any extra expense.

CLAIMS ASSISTANCE - Call Picarro, Inc. Customer Service or your local distributor for assistance. Give our representative the full details of the problem. Helpful information or shipping instructions will be provided. If requested, estimates of the charges for non-warranty or other service work will be supplied before work begins.

RETURN PROCEDURE - Customers must obtain a Return Merchandise Authorization Number from Picarro, Inc. prior to returning units. Products being returned for repair must be shipped in their original shipping cartons to avoid damage.

APPENDIX A – REMOTE DATA ACCESS

A.1 Picarro Serial Communication

The analyzer supports an RS-232 physical command interface, which can be used to control the instrument and to retrieve concentration data. Not all features of the instrument are available on the serial interface. For details on how to use the serial command interface, please see the Programming Guide (included in pdf format on the installation CD). This command set may also be used across a TCP/IP interface through an Ethernet connection. Please contact Picarro for further details.

A.2 Remote Data Access

Using the RemoteAccess.ini file, the analyzer can be configured to automatically:

- Send data from the instrument to a list of e-mail accounts.
- Measure the offset of the host computer system clock from a set of Internet time servers and (optionally) to resynchronize the clock based on this information.

The Internet connection need not be permanent and may be a dial-up connection accessible via a user-supplied USB modem. The task of sending data and/or synchronizing the clock on the analyzer is performed using the C:\Picarro\G2000\HostExe\RemoteAccess.exe program. This program can be set up to run periodically using the Windows task scheduler at a user-configurable frequency. If a dial-up connection to the Internet is employed, it is used only on demand to minimize the connection time.

Each time that the RemoteAccess.exe program runs, it appends information to a log file, which keeps a record of the results of the time synchronization and of the files sent by e-mail. The RemoteAccess.exe program is configurable by means of an initialization file, which includes information such as the login credentials for the dial-up connection, the e-mail account and the list of time servers.

The initialization file is:

C:\Picarro\G2000\AppConfig\Config\UtilitiesRemoteAccess\RemoteAccess.ini

and it should be placed in the same directory as the executable RemoteAccess.exe. The file has one required section named LOGGING and three optional sections named NTP, DIALUP, and EMAIL. The logging section has a single key Logfile whose value is the path to the log file. Once this log file exceeds 64K bytes in length, it is backed up, appending a numeric extension to the file name, and a new file is opened. A total of ten backup log files are kept.

NTP

The NTP section controls querying the Internet time servers using the SNTP protocol (RFC4330) and the resetting of the clock on the host computer. If the section is not present, time synchronization is not carried out. The keys Server1, Server2, etc., are used to specify the URLs of the time servers. If the UpdateClock key is set to "true," the offset is applied to the host clock. Otherwise, the offset is recorded, but the host clock is not changed.

Email

The EMAIL section controls the sending of the data files as e-mail attachments. If the section is not present, e-mail messages are not sent. The key Directory specifies the directory that contains the data files. When the program is run, files in this directory are sent to the specified recipients and the files are deleted. To avoid problems with incomplete files, programs that place files into this directory should do so using an atomic operation, such as a rename. The Server key is set to the name of an RFC2821-compliant SMTP server that sends the e-mail messages.

The FROM key is the e-mail address from which the messages are sent. Note that some SMTP servers check that the source is permitted to send email while others allow any name in this field. The collection of e-mail addresses to which copies of the e-mail is sent is specified by the keys To1, To2, etc. The Subject key is used to fill the subject field in the email header and may be set to any string. Depending on the SMTP server, it may be necessary to use authentication before e-mails can be sent, as described in RFC2554. If such authentication is not needed, the key UseAuthentication is set to false. If this key is set to true, two additional keys Username and Password must also be specified for the e-mail account.

Dial-up

The DIALUP section is used if a dial-up connection to the Internet needs to be established when the program runs. If the section does not exist, a permanent connection is assumed to be available for carrying out the other tasks specified in the initialization file. The connection name key specifies the name of the dial-up connection to use, as listed under Network Connections in the Control Panel. The values of the keys Username, Password and Number are used to make the connection.

Examples of “RemoteAccess.ini” File

[LOGGING]	Logfile=c:/temp/RemoteAccessLog
[NTP]	Server1=time-a.nist.gov Server2=time-b.nist.gov Server3=time-a.timefreq.bldrdoc.gov Server4=time-b.timefreq.bldrdoc.gov Server5=time-c.timefreq.bldrdoc.gov Server6=time.nist.gov Server7=time-nw.nist.gov UpdateClock=1
[DIALUP]	ConnectionName=Picarro Dialup Access Username=user Password=password Number=14085551212
[EMAIL]	Server=smtp.servername.org Directory=c:/picarro/mailbox From=instrument@picarro.com To1=recipient1@site1.com To2=recipient2@site2.com Subject=CRDS data from Silverstone instrument UseAuthentication=0

APPENDIX B – THE DATA FILE VIEWER

B.1 Introduction

The Picarro Data File Viewer software is located in the Picarro Utilities folder on the desktop. The software allows the user to graph and to conduct statistical analysis of the raw data. Additional functions include Allan Variance plot and quadratic or polynomial fittings.

The Picarro Data File Viewer includes two main menus: **File** and **New**.

Figure 50:
File Menu


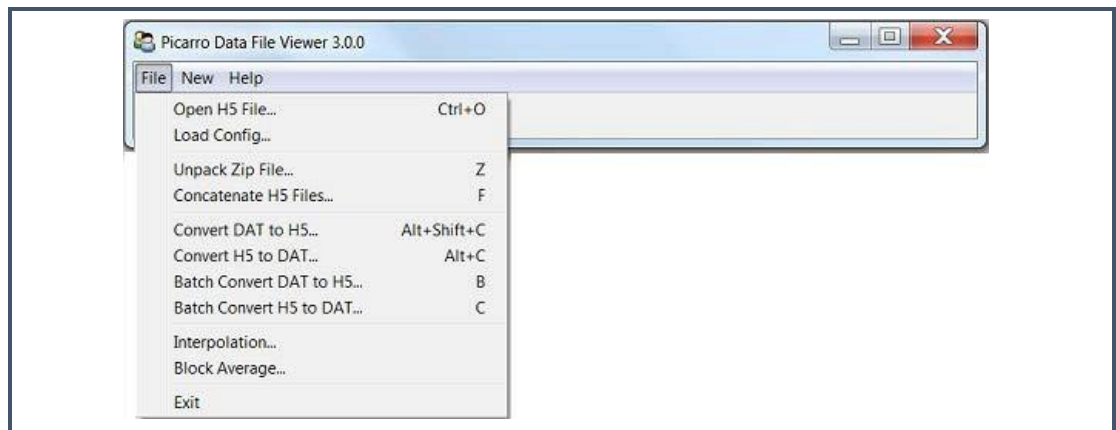


Figure 51:
New Menu

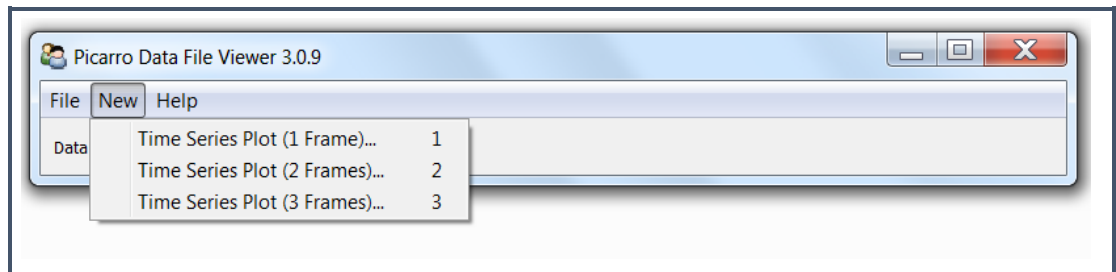

B.2 The File Menu

This section describes the functions available from the Data File Viewer File menu.

Open H5

File > Open H5 opens a Picarro data file (HDF5 format) for data analysis and visualization. After opening the data file, you can create a new time series plot. Refer to *New Time Series Plot* for more information.

Load Config

File > Load Config loads a configuration file (ini format) to restore parameters of a workplace. Refer to *Save Configuration* for more information.

Unpack Zip File


Use **File > Unpack Zip File** to concatenate all H5 files inside the zip file into a single H5 file. Refer to *Concatenate H5 Files* for details.

Concatenate H5 Files

Use **File > Concatenate H5 files** to concatenate multiple files and zip archives of H5 files into a single H5 file. Navigate to the desired folder, or use the **Define date range** button to specify a date range of files to concatenate. (See next section.)

After selecting the path of the data files, Data File Viewer will automatically search an H5 file in the specified zip/folder and look for all available variables in the H5 file. The variables are then listed in the left panel, and users can use “>>” button to move variables to the right panel for concatenation.

Figure 52:
Select Variables
Form

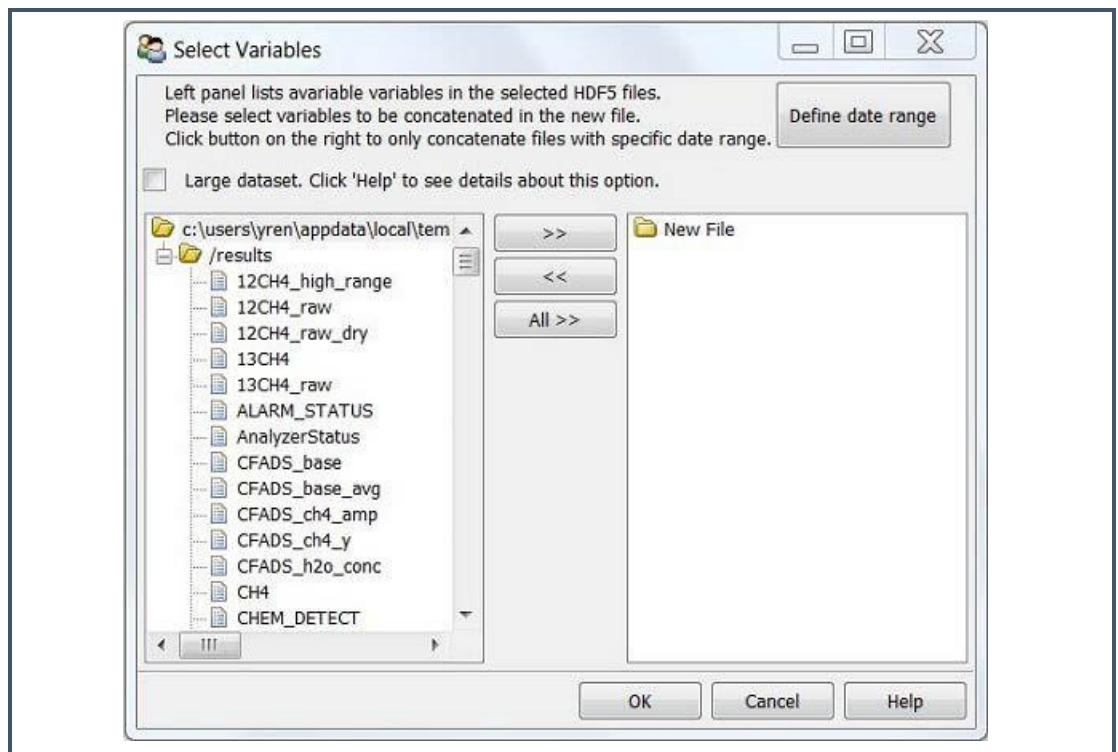

Define Date Range

Data File Viewer can search data files within the desired date range and then concatenate such files into an H5 file.

By default, TimeZone is set to your local time zone. However, if data were taken elsewhere, select the time zone where data were taken.

Select File > Concatenate H5 files, then click Define date range to specify the desired date range.


Figure 53:
Defining Date Range



Data files are saved in directory trees named by date and time option.

Picarro software saves data in a directory tree that is named by the creation year, month, and day. (See example directory tree in the following image.) Select this option if the target folder has this file structure. This way, Data File Viewer will only search folders within the desired date range, which can substantially reduce processing time.

Figure 54:
File Structure of Data File Viewer

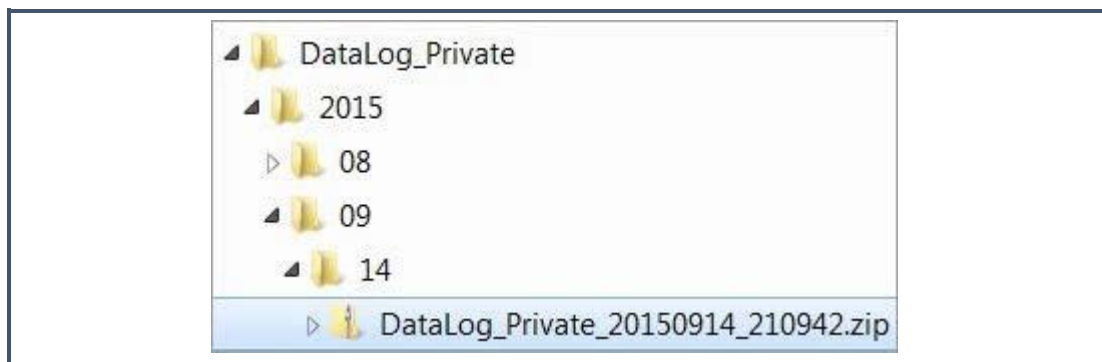



NOTE

To save processing time, Data File Viewer does not open data files, but only determines data acquisition time based on the file name.



CAUTION

Do not define a time range for data files whose names have been changed.



NOTE

Data File Viewer does not concatenate data files exactly within the defined time range. This is because the time extracted from file name is different from the data acquisition time. To not miss data points, Data File Viewer expands the specified time range, so the resulting dataset normally has a wider time range than the user specification.

Convert DAT to H5

Select **File > Convert DAT to H5** to convert a file in DAT format to HDF5 format. These formats are described below:

- DAT format: DAT files accepted by DatViewer store tabular data (numbers and text) in plain text.
- Each line of the file is a data record. Each record consists of one or more fields separated by whitespaces.
- The first line of data file indicates column names.
- There must be a field “EPOCH_TIME” to store the acquisition epoch time (expressed as seconds since Jan 1, 1970) of the data. Otherwise, the first and second fields must be “DATE” and “TIME”. The “DATE” field must have the format “mm/dd/yyyy” or “yyyy-mm-dd”, and the “TIME” field must have the format “HH:MM: SS(.sss)” where (.sss) is an optional fraction of seconds.
- HDF5 format: HDF5 is a data model, library, and file format for storing and managing data. (See the HDF5 Home Page on the HDF Group Web site for more information.) When converting DAT to HDF5 format, Data File Viewer creates a table named “results” to the contained data.

Convert H5 to DAT

Select **File > Convert H5 to DAT** to convert a file in a HDF5 format to DAT. These formats are described in *Convert DAT to H5*.



NOTE

Data File Viewer does not concatenate data files exactly within the defined time range. This is because the time extracted from file name is different from the data acquisition time. To not miss data points, Data File Viewer expands the specified time, column name “fineLaser-Current_1_controlOn” will be replaced with “fineLaserCurr_1_ctrlOn”.

Interpolation

Interpolation describes the method for constructing data points with a range of a discrete set of known data points. Select **File > Interpolation** to perform interpolation on a time grid with a constant interval.

Block Average

Select **File > Block Average** to divide a dataset into small blocks based on a user-defined block size. The average is calculated for data in each block, and the results are saved in a new H5 file.



NOTE


The specified block size must be greater than the average data interval.

Because the data interval is normally not a constant (unless interpolation is performed), fluctuations in the data interval will affect block averaging if the block size is comparable to the average data interval.

B.3 New – Time Series Plot

You can specify to include create time-series plots with one, two, or three frames. New plots display in the Time Series Viewer.

Figure 55:
Time Series Viewer

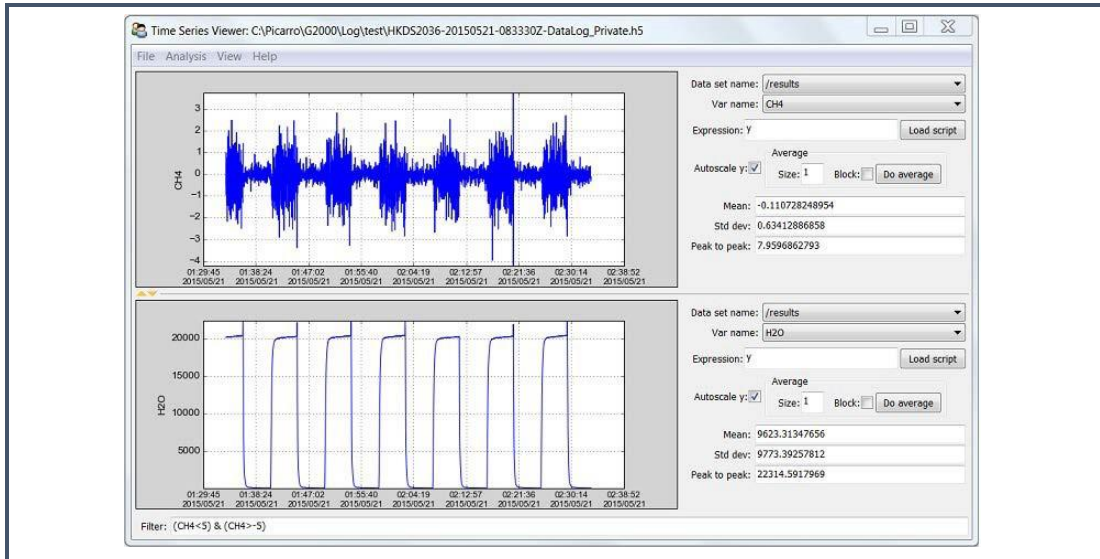


The next section describes the options available on the Time Series Viewer menu bar. Refer to *Time Series Viewer Canvas* on page 87 for information the Time Series Viewer UI features and options.

B.4 Time Series Viewer Menus

The Time Series Viewer form includes the following menus:


Figure 56:
Time Series Viewer Menus

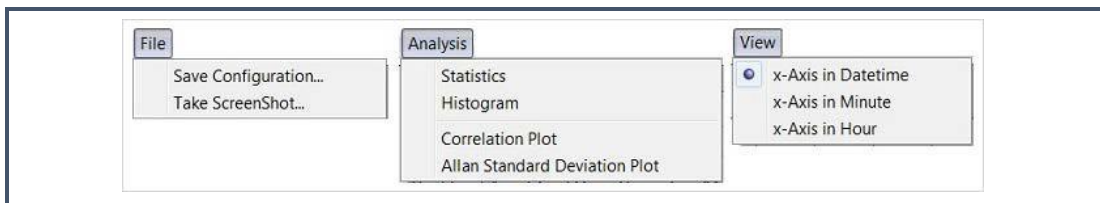

Time Series Viewer File Menu

Use the File menu to save a configuration or take a screenshot.


Figure 57:
Time Series Viewer – File Menu

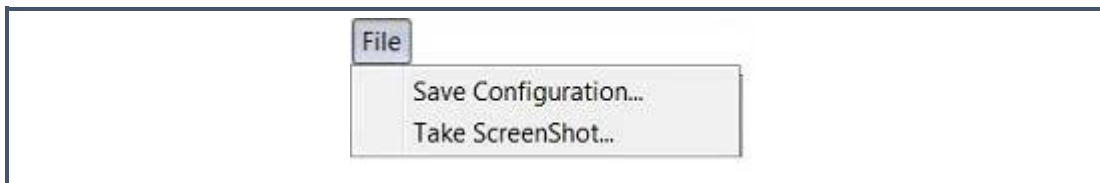

Save Configuration

Click **File > Save Configuration** to open the Feature Capture form. With this form, you can save figure properties, expressions, filters, and other settings to a configuration file so that it can be easily loaded in the future.

Figure 58:
Time Series
Viewer –
Feature Capture

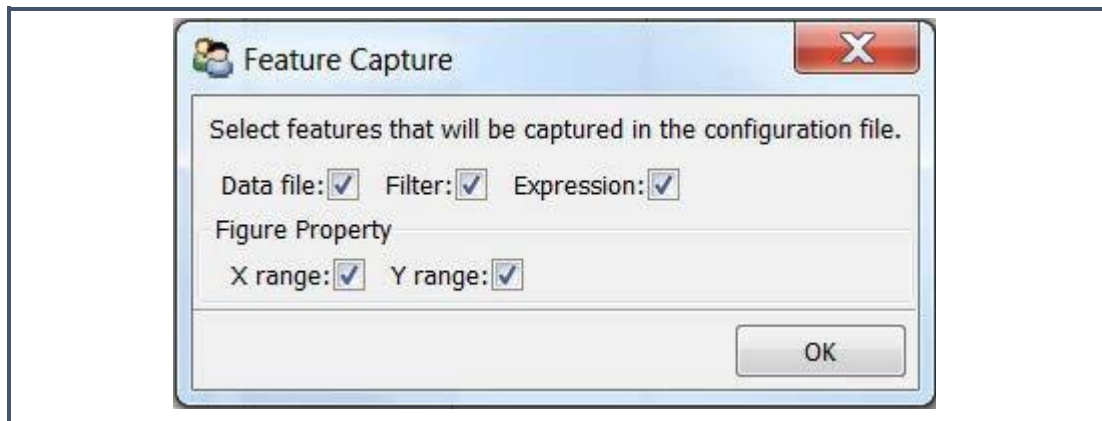

CAUTION

If a feature is not captured, it will be omitted when the configuration file is loaded.

Depending on the features captured, loading a configuration file can have different effects. For example:

- If all features are captured, a saved workplace is reproduced.
- If Data file is not captured, saved parameters will be applied to the data file in memory.
- If Expression is not captured, plots will not be transformed.
- If X (Y) range is not captured, figures will be auto-scaled on the x (y) axis.


Take ScreenShot

Use **File > Take ScreenShot** to take a screenshot of the Time Series Viewer and save it as a .png to a specified file.

Time Series Viewer Analysis Menu

Use the Analysis menu to calculate statistics, generate a histogram, and to plot correlations and Allan Standard deviations.

Figure 59:
Time Series
Viewer –
Analysis Menu

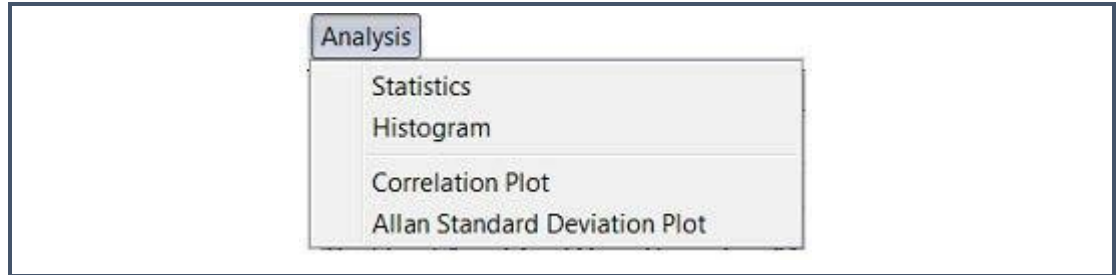


Statistics

Use **Analysis > Statistics** to calculate mean, standard deviation and peak to peak for all plots in the current window.

Histogram

Use **Analysis > Histogram** to generate a histogram of data as shown below.

Figure 60:
Histogram
Window

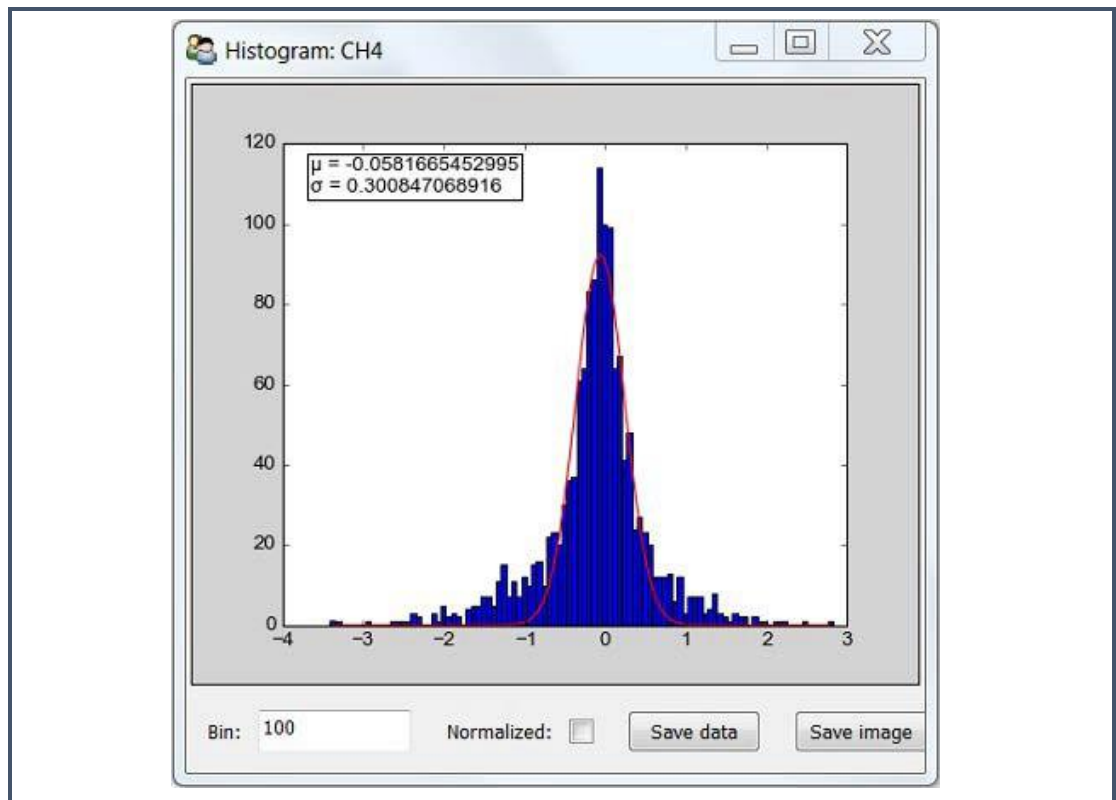

- **Red Line:** A Gaussian function fitted to the histogram. Fitting results of μ and σ are shown in the top-left corner of the plot.
- **Bin:** Specifies the number of intervals that the range of values is divided into.

- **Normalized:** When selected, the sum of the histograms is normalized to 1.
- **Save data:** Saves histogram data to a CSV file.
- **Save image:** Saves the histogram image as a JPEG/PNG/PDF file.

Correlation Plot

Use **Analysis > Correlation Plot** to plot Y-axis data in one frame versus that in the other. This can be used when two or more frames exist in the current Time Series Plot window. See Correlation/XY Plot on Page 89 for details.

Allan Standard Deviation Plot

Use **Analysis > Allan Standard Deviation Plot** to create an Allan Standard Deviation plot (versus a standard deviation plot) for data in the current window. See [Allan Variance](#) for more information.

Time Series Viewer View Menu

Use the View menu to view X-axis information in date-time, minute, or hour format.

Figure 61:
Time Series
Viewer –
View Menu




NOTE

When switching from DateTime to Minute or Hour, the X-axis data is subtracted from the earliest point shown in the panel and then converted to the desired unit.

The Time Series Viewer Canvas

The Time Series Viewer canvas is comprised of interactive graphs and a variety of configuration options.

Figure 62:
Time Series
Viewer Canvas

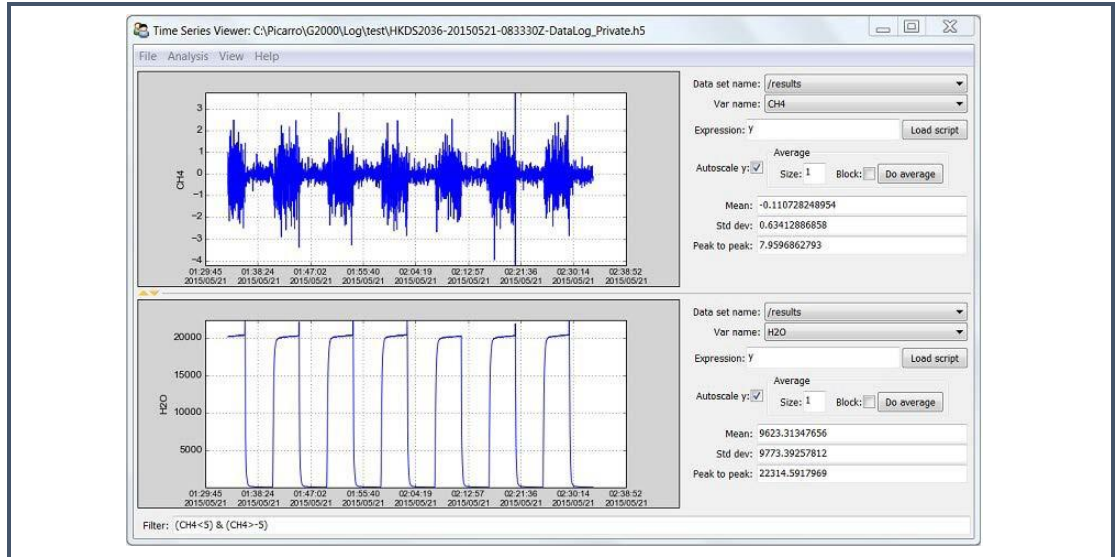

Mouse Options and Graph Transform


The following mouse actions can be used in the canvas graphs:

- Left click and drag: Zooms into the selected area of the plot.
- Left click and drag with the SHIFT key down: Pans the plot.
- Left click and drag with CTRL key down: Zooms out from the plot.
- Left click and drag with ALT key down: Stretches the plot.
- Right-click: Opens an additional menu. Refer to the next section.

Right-click Menu

Right-clicking on the canvas provides opens a pop-up menu.

Figure 63:
Canvas Right-
click Menu

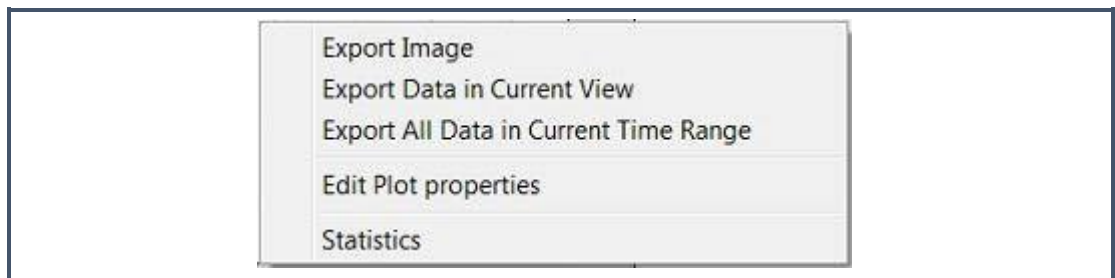

Export Image: Exports the current plot as a jpeg, png, or pdf file.


Export Data in Current View: Exports only date/time and the selected variable in the current view to an HDF5 or CSV file.

Export All Data in Current Time Range: Exports all variable columns of the selected dataset in the current time range to an HDF5 file. Refer to *Concatenate H5 Files* on page 79 for more information.

Edit Plot properties: Opens the Image Editor form, where the following options can be specified:

- **Title:** Edits the title of the plot.
- **Line:** Specifies the line pattern of the plot. If **None** is selected, the data points will be plotted without connecting lines.
- **Marker:** Specifies the marker type to indicate data points. If **None** is selected, data points will not be shown.
- **Min and Max:** Specifies the minimum and maximum of data range for the X-axis and the Y-axis.
- **x[0]:** Sets the earliest time of the dataset as the minimum of the X-axis.
- **Time zone:** Sets the time zone for date/time variables. This defaults to the local time zone.
- **Label:** Specify a label for the X-axis and the Y-axis.

Figure 64:
Image Editor
Form

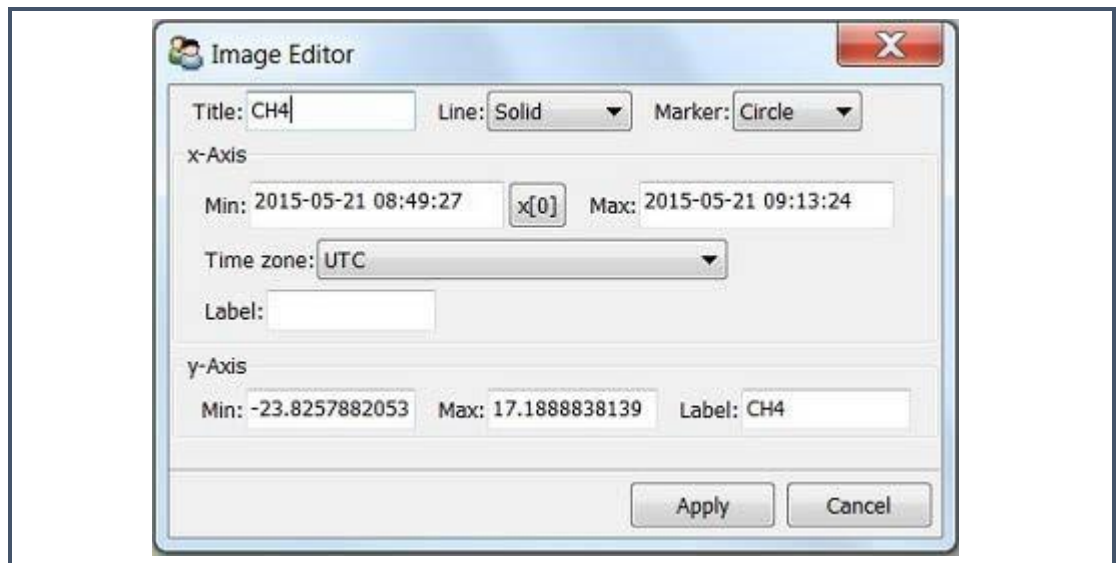

- **Statistics:** Calculates the mean, standard deviation, and peak to peak for data in the current view.

Dataset Name and Var Name

An HDF5 file can store one or more tables. Each of these tables is called a Dataset. A table can contain one or more columns. Each column is called a variable (Var).

Use the **Dataset name** drop down to select the dataset that will be used for this time series graph. Use the **Var name** drop down to select the column in the dataset to use in the graph.

Autoscale Y

When the Autoscale Y option is selected, the Time Series Viewer will autoscale on the Y-axis to make sure that all data within the range of the X-axis is displayed.

Average

If **Block** is selected, a block average is calculated when you click the Do average button is clicked. Otherwise, a moving average is calculated.

For a block average, **Size** specifies block size in unit of a minute. For a moving average, **Size** specifies subset size in unit of data points.



REMINDER

Averaging is performed after the Filter and Expression are performed.

Mean, Std Dev, and Peak to Peak

The **Mean**, **Std dev** (Standard deviation) and **Peak to peak** fields provide all the statistical information of data in the current view.

Correlation/XY Plot


The Correlation/XY Plot window includes two menu items: **File** and **Analysis**. For details about the File menu, see *Save Configuration* on page 84.



REMINDER

The canvas in this plot is interactive. For details about the plot canvas, see *The Time Series Viewer Canvas* on Page 87.

Figure 65:
Correlation/XY
Plot

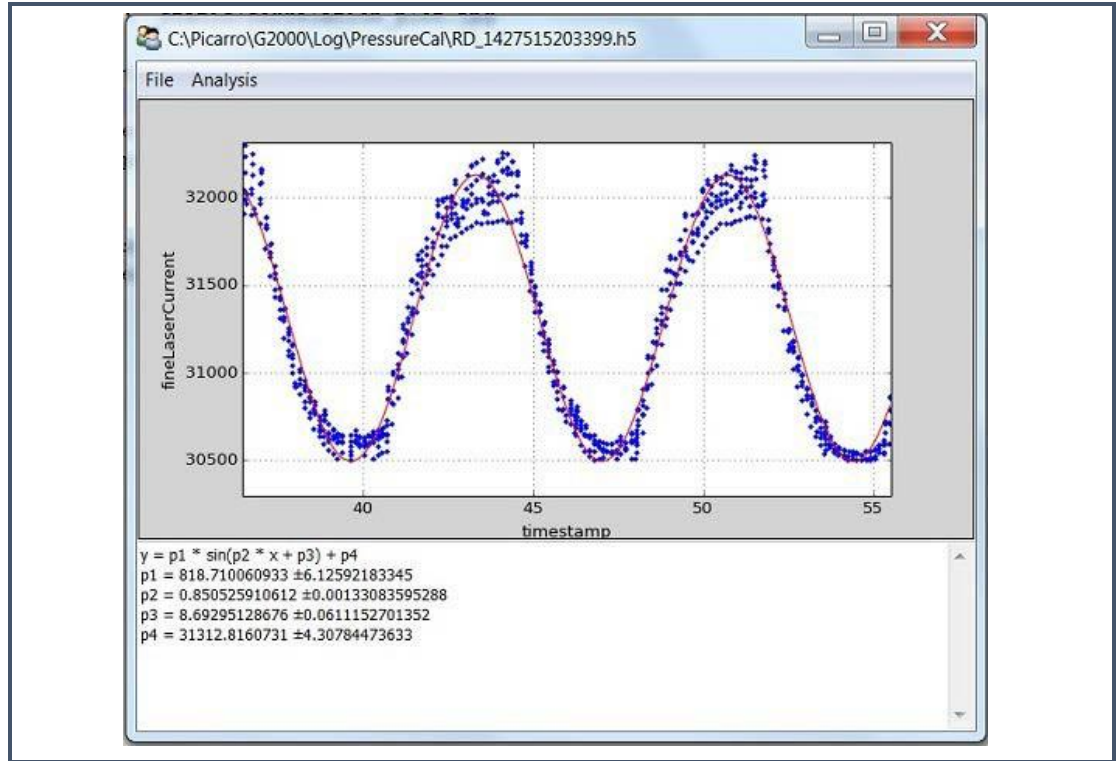


Analysis Menu

The Analysis Menu includes three options: Fitting, Integration, and Statistics.

- **Fitting** allows you to specify one of four fitting methods to include in the Correlation/XY plot:
 - i. **Linear fit:** Specifies to fit to linear function $y = c_1x + c_0$
 - ii. **Quadratic fit:** Specifies to fit to quadratic function $y = c_2x^2 + c_1x + c_0$
 - iii. **Polynomial fit:** Specifies to fit polynomial function of degree n : $y = \sum c_n x^n$
 - iv. **Curve fit:** Specifies to use non-linear least squares to fit an arbitrary function to data.


Figure 66:
Analysis
Menu

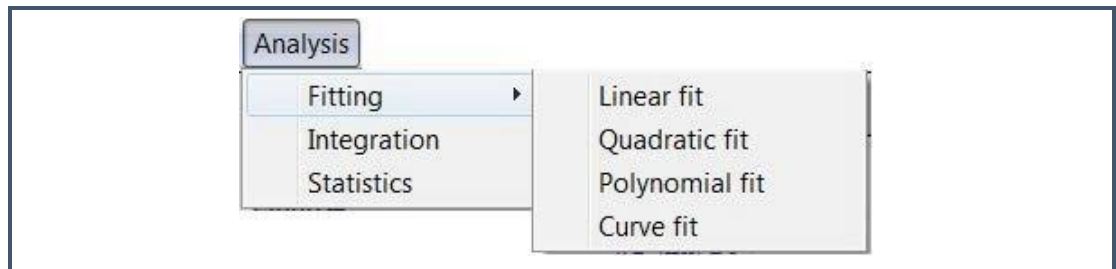

- **Integration** calculates area under the curve using the composite trapezoidal rule.
- **Statistics** calculates mean, standard deviation, and peak to peak for data in the current view.

After applying any of the above Analysis options, the results, statistics, or fitting function with coefficients are displayed in the lower portion of the Correlation Plot window.

Figure 67:
Results of
Quadratic
Fitting

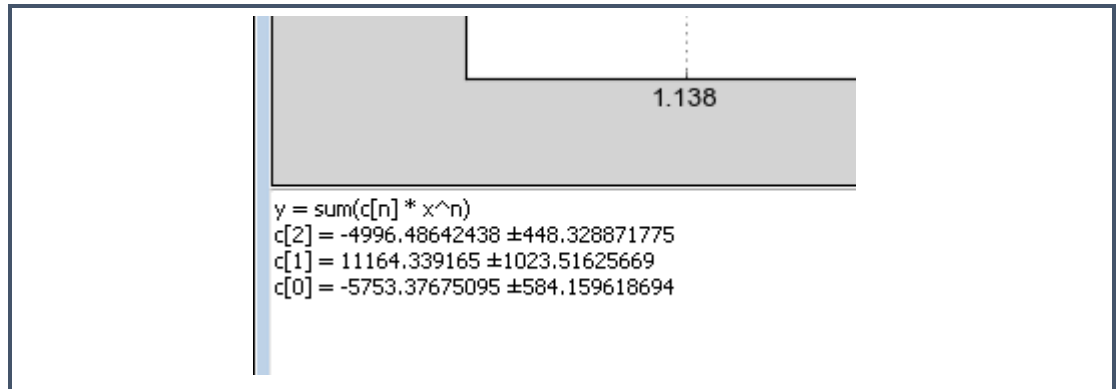

APPENDIX C – EXTERNAL VALVE SEQUENCER

C.1 Introduction

The Picarro analyzer can control two types of valves:

- **Solenoid valve(s):** DC voltage powered valve with normally open (NO) and normally closed (NC) positions. These can be either 2-way or 3-way valves.
- **Rotary Selector valve:** digitally controlled valve used to send selected flow from one of many inputs (up to 32) into the analyzer

Both types of valves can be simultaneously controlled through a common software interface called the 'External Valve Sequencer,' which is available from the Tools menu in the GUI.

C.2 Default Configuration

For all models of Picarro analyzers, the rotary valve control is disabled in the factory default setting. The solenoid valve control, however, is ready to use by default for all solenoid valve connectors.

C.3 Setting Up Solenoid Valves

The Valve Sequencer software can control up to six solenoid valves. Each valve should operate using 12 VDC with a current requirement of <1.5 Ampere maximum. This analyzer comes with a cable that can be connected to the solenoid valves.


The valve connector cable should be connected to the 15-pin connector at the lower left corner of the analyzer. There are six pairs of wires with connectors labeled V1, V2, ... V6 with 2-pin female Molex connectors (Molex#43020-0200) for connection to the solenoid valves. For valves wired with matching Molex connectors, connect V1 to the solenoid valve 1, V2 to solenoid valve 2, etc. Do not connect the solenoid valve to the analyzer ground -- use only the provided electrical connectors.

C.4 Setting Up Rotary Selector Valve

A multi-position rotary selector valve can be controlled by the Valve Sequencer software. It is controlled by standard serial commands in the Valco (VICI) protocol. Valco rotary valve models SD, SC, SF, ST, and STF are all supported. However, not all configurations will be appropriate, due to tube diameter, pressure or materials of construction. This setup will also require a Valco micro electric high torque actuator. A single combination package such as EMT2ST16MWE includes a 16 position, low- pressure ST valve in stainless steel, 1/8" tubing, 2" standoff, and micro electric high torque actuator. Please refer to http://www.vici.com/vval/st_8-1.php#16pos for further options and consult with Valco directly for more details.

The Valco controller should be installed per manufacturer's instructions. The 9-pin female, connector cable (female) should mate with its corresponding, male port of the analyzer and labeled "MPV." Please note the 9-pin connector cable is not supplied with the instrument.


Figure 68:
8-port Valve Setup

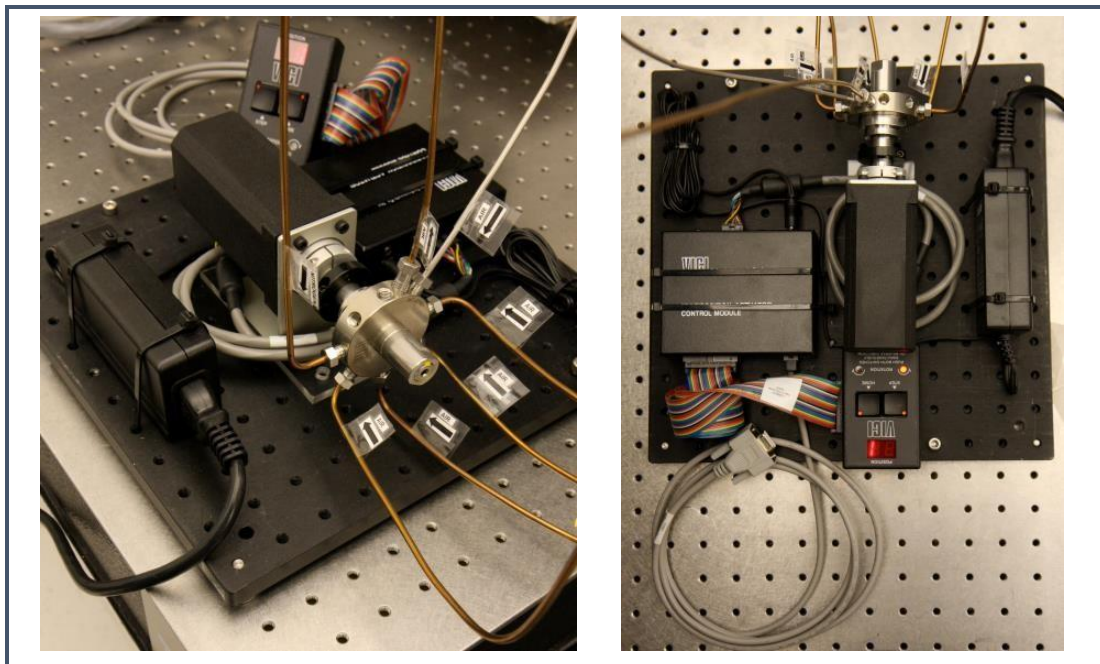

C.5 Valve Sequencer Software

The Valve Sequencer software allows the user to set steps in which solenoid valves are turned on/off and the rotary selector valve is set to a single position. The Picarro valve sequencer window appears below:

Figure 69:
Valve
Sequencer GUI

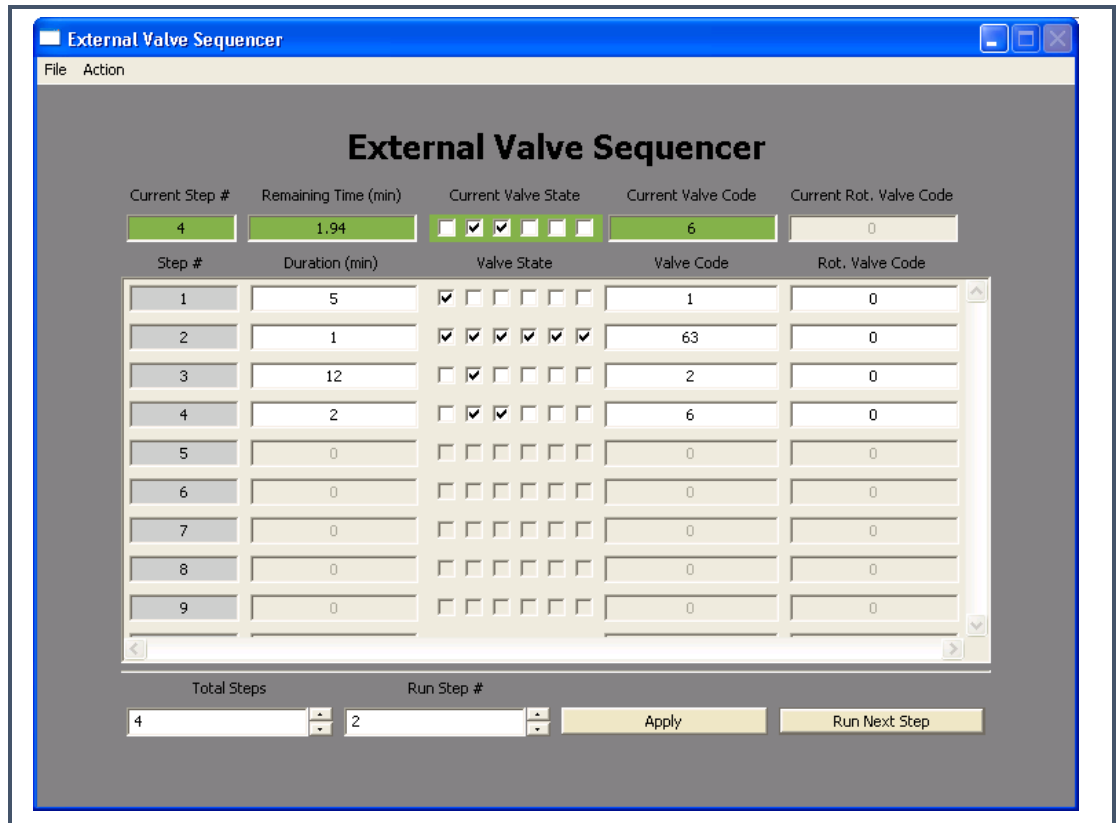

The current step, elapsed duration, and valve states are shown in the topmost row of the valve sequencer command window. The duration of each step is set in minutes; for example, 1 minute and 30 seconds correspond to 1.5 minutes.

Please note the number of steps correlates with the total steps in the sequence, and the count of the steps starts at “0.” The first step in the sequence is designated Current Step “0,” and the second step in the sequence is designated Current Step “1.” The “Go to First Step” menu item under “Action” restarts the sequence from step 1.

Different sequences can be created and saved in the software. Use the “Sequence File #” field to select which file the sequence will be saved to (0 to 10 are the available choices). Click the up/down arrows to select the desired number.

C.6 Configuring a Valve Sequence

Each “step” sets the rotary valve to a single position and activates the indicated solenoid valve(s) for a set period. Multiple steps can be carried out in sequential order to switch between different gas sources, flush out a manifold, or other gas handling operation. Create the number of desired steps in the sequence by clicking the up/down arrow for “steps.”

For each step, select the box for each solenoid valve to be opened. The checkmark in the “current valve state” window indicates a solenoid valve is open. Note: In this example, we assume normally closed (NC) valves are used. A check indicates current is flowing to the valve thus powering it open. The positions from left to right correspond to solenoid valves V1...V6.

The rotary selector valve position can be set in the column labeled “Rot. Valve Code.” Enter the number that corresponds to the desired valve position. A value of 1 in this field corresponds to position 2 on the Valco valve. Only one rotary position can be selected per step.

Step duration is determined by the value entered in the “duration (min)” field, where the duration of the step is in minutes. If duration values are set to <0.1 minutes, they may not be carried out correctly.

The “valve code” field is a configuration- dependent, read-only display field that shows the total state of that particular step in a numerical code. Should the most upper right gray box display a value of 512,256 or be grayed-out, either no rotary selector valve is connected to the instrument, or the valve is not functioning. For each individual measurement, the analyzer makes, the valve codes and rotary valve positions corresponding to the valve state(s) at that point in time are saved alongside the concentration data.

Once the valve sequence has been programmed, it can be saved using the button “Save Valve Sequence File.” The sequence will be saved under the sequence file number selected.

C.7 Loading and Running a Saved Sequence

To load a valve sequence file, select the desired “sequence file number” and press “Load Valve Sequence File.” If the user has been running a different sequence from the one that was loaded, the user needs to press “Next Step” to initialize the newly selected sequence.

To run a sequence file, press “Enable Sequencer.” This button will turn to “Disable Sequencer” once the sequence starts. (The sequencer should be

activated if it was disabled, but not necessarily to change from one sequence to another.) The sequence will repeat itself indefinitely until disabled or the software is exited. If enabled, the sequence will continue to run after the “close sequencer window” button is pressed.

If desired, the valve sequence can be forwarded to the next step of the sequence by pressing the “run next step” button. To stop the sequencer file, use the “Start/Stop Sequencer” menu item under the “Action” menu. This will leave all valves in their current state. In some situations, it is convenient to program the last step in the sequence to be a safe or default valve state. The sequencer can be advanced to the last step should the user need to put the solenoid or rotary valves into a safe/default state. The “Reset All Valves” de-activates all valves. Using the “Hide Sequencer Interface” closes the window, but if the sequencer is enabled, it will continue to run in the background. To jump to a particular step, increment the “run step 3” and click “Apply.”

Both solenoid and rotary valve codes are recorded in columns in the output data files indicating the active valve configuration respective to when data is taken. These codes can be used as event timing flags. For example, if no solenoid valves are present, the codes will be recorded regardless of whether a valve is connected or not.